

## Long-slit Spectroscopic Data Reduction For 90-inch Bok Using The Image Reduction Analysis Facility

SARA E. HOLEMAN<sup>1</sup> AND SANCHAYEETA BORTHAKUR<sup>1</sup>

<sup>1</sup> *School of Earth & Space Exploration  
Arizona State University  
781 Terrace Mall, Tempe, AZ 85287, USA.*

### 1. INTRODUCTION

This manual is specific to the procedure of reducing galaxy long-slit spectroscopic data from the Boller & Chivens Spectrograph aboard the 90" Bok telescope at Kitt Peak National Observatory. The Image Reduction Analysis Facility (IRAF) is utilized to complete the steps of image reduction. This manual will discuss methods of removing instrumental signatures, correcting for radiation events from environmental conditions, and the data frames for wavelength calibration for the particular data presented. We refer the reader to *A User's Guide to CCD Reductions with IRAF* [1] for more general and in-depth information.


In Section 2, the spectra data format and calibration frames will be introduced. Section 3, will present IRAF tasks with their appropriate parameters following the steps of the initial reduction. Lastly, Section 4 will highlight the importance of wavelength calibration with the tasks that create a two-dimensional function to project the scientific data on a linear wavelength scale.

### 2. SPECTRA FORMAT & CALIBRATION FRAMES

Upon collecting data, several artificial components are superimposed on the CCD due to instrumental factors such as the electronic bias level, the multiplicative gain and illumination variation, and any non-linear radiation events, particularly cosmic rays. Therefore, removing these fallacious data components is required before any scientific data analysis.

#### 2.1. Two-dimensional galaxy spectra

Analyzing a two-dimensional spectrum, you will see the flux as a function of wavelength. Reference Figure 1 for the slit position relative to the major axis of the galaxy and reduced two-dimensional spectra of a low-redshift galaxy, showing prominent H $\alpha$  emission lines extending out to the galaxy disk.

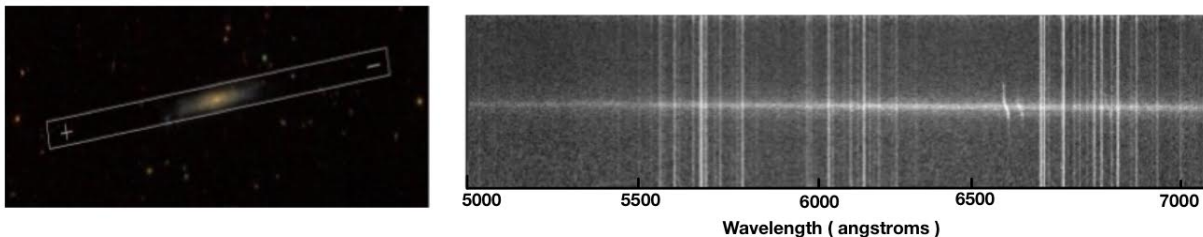

**Figure 1.** Left: Panel showing the position of the slit with respect to the target galaxy. Right: The two-dimensional data showing the galaxy continuum. The H $\alpha$  emission line can be clearly identified thus indicating the rotation of the galaxy by its red and blue shifted wavelengths. The vertical lines along the frame are emission skylines.

#### 2.2. Overscan Region & Bias Frames

While obtaining the data, an additional voltage is added to the output signal, that comes from the CCD's sensor while reading out the signals from light exposure. Therefore, the additional counts need to be subtracted from all the frames as the initial step of the image reduction. This step is referred to as "bias subtraction" and leads to identifying the overscan region where no incident light falls onto it during exposures. The overscan region, which exists at the tail end of the CCD (around 30 pixels in length and the full width of the CCD), allows only the biased voltage to be defined. As this value varies slightly along the columns of the CCD, a first-degree polynomial fit is easy to obtain or simply averaging the data over the columns of the overscan region is sufficient as well. The overscan region will be trimmed off after subtracting the fit from your images.

It is typical to expose the CCD to a flash of uniformed light (pre-flash) to ensure no charge loss during signal transfer. This superimposes yet another number of counts on the data that the bias subtraction may not account for. This can be corrected by combining a collection of bias frames. These frames are "zero" second integration exposures, and by median combining to a single frame, the small pre-flash counts can be subtracted from the other reduction frames.

### 2.3. Dark Current

Dark current is a temperature-dependent source of additional noise accumulating over long exposures. In the case of the data presented here, the dark current is not significant and, therefore, not corrected. However, it is possible to check if correction is needed. Because the overscan and bias levels would have already been removed from all calibration frames, if any counts are visible on the dark frames, then that will be a significant source (i.e. real) of dark current.

### 2.4. Flat field correction & Fixing non-linear pixel sensitivities

Flat field data is fundamental to help define light sensitivities along the CCD, it is an image of uniformly distributed light. Dividing the raw data by a normalized flat field, (done in an IRAF), corrects for uneven illumination, dust, and vignetting in an image. The preservation of counts is important for spectra analysis. Thus, two additional requirements must be met before normalization. Firstly, it is useful to interpolate over any bad pixels on the CCD, which ultimately appear in the images. Interpolation can be done by creating a bad pixel mask that will be used during the final correction of the reduction process. Second, possible wavy distortions may appear on the images showing a large-scale variation in light that is wavelength-dependent caused by several potential factors like the temperature or color imbalances of the flat-field lamp. Two tasks will be used for both processes and introduced in the following section.

## 3. IMAGE REDUCTION WITH IRAF

IRAF is a command language program with packages that perform tasks depending on the specific parameters given. This section will cover image reduction package tasks and instructions on the parameters that are used to correct raw data. To edit the parameters of a task enter **epar** in the command line followed by the task name you'd like to edit (e.g., **epar zerocombine**). Once the parameters are appropriate to the step at hand, save them by **control-Z** on your keyboard.

### 3.1. Characterizing the data

To familiarize yourself with the data, use the task **cdlist** to analyze the information of individual frames present in the folder you are working in. The output can be seen in Figure 2, where each frame has header information associated with the object type defined by the observer at the telescope. A translation table regarding the header information must be set up so that IRAF can reference the appropriate data when running called tasks. The **setinstrument** task will do so and is located in the **ccdred** package. For spectral reduction, the input parameter needs to be set **specphot**. Upon running this task, you will enter another parameter editor for the **ccdred** package. Note that the primary defaults of the package will work just fine but ensure that the pixel type parameter is set to **real**.

```


ccdred> ccdlist
CCD images to listed (*.fits):
b0001.fits[1220,326][ushort][flat][]:flat
b0002.fits[1220,326][ushort][flat][]:flat
b0003.fits[1220,326][ushort][flat][]:flat
b0004.fits[1220,326][ushort][flat][]:flat
b0005.fits[1220,326][ushort][flat][]:flat
b0006.fits[1220,326][ushort][flat][]:flat
b0007.fits[1220,326][ushort][flat][]:flat
b0008.fits[1220,326][ushort][flat][]:flat
b0009.fits[1220,326][ushort][flat][]:flat
b0010.fits[1220,326][ushort][flat][]:flat
b0011.fits[1220,326][ushort][flat][]:flat
b0012.fits[1220,326][ushort][flat][]:flat
b0013.fits[1220,326][ushort][flat][]:flat
b0014.fits[1220,326][ushort][flat][]:flat
b0015.fits[1220,326][ushort][comp][]:arc_J1120
b0016.fits[1220,326][ushort][object][]:J1120
b0017.fits[1220,326][ushort][object][]:J1120
b0018.fits[1220,326][ushort][object][]:J1120
b0019.fits[1220,326][ushort][object][]:J1120
b0020.fits[1220,326][ushort][comp][]:arc_J1120
b0021.fits[1220,326][ushort][comp][]:arc_J1122

```

**Figure 2.** `ccdlist` will output a list shown above presenting the header information for each frame in the dataset regarding its dimensions with its data and object type.

### 3.2. Locating the overscan region


Examining the flat field data can be useful in identifying the location of the overscan region. There are two ways of defining this region. First, by visual inspection and locating the specific column where the cutoff is along the CCD and, secondly, by examining the plot of the flat field with the IRAF `implot` task. Figure 3 presents the two forms we can visualize the data.

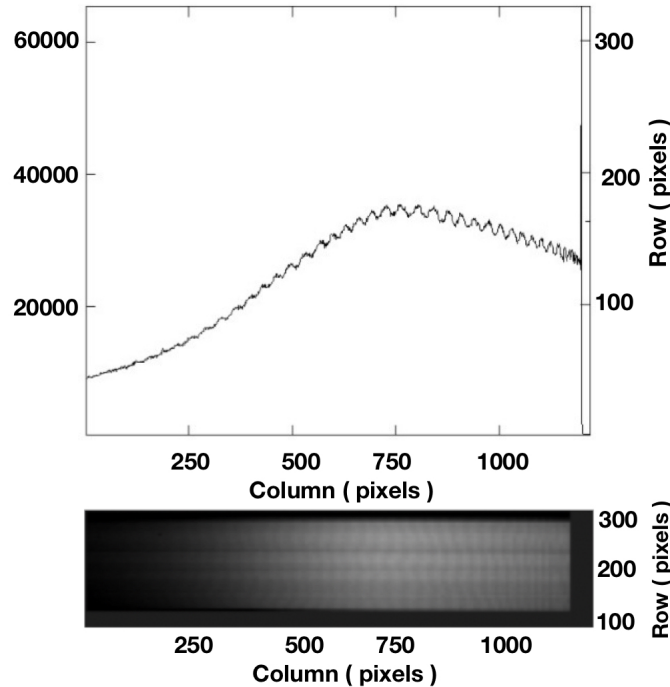

**Figure 3.** Top: Flat field plot showing the light variation across the CCDs pixel columns. Bottom: For visual inspection, the overscan region is the dark portion of the CCD where there's no penetrating light (see far right, top, and bottom of the frame).

### 3.3. Combining Bias Frames with *zerocombine*

To ensure no charge loss during signal transfer, an additional number of counts are added to the data by pre-flashing the CCD before exposures. Combining data frames is a technique that is frequently used in this manual. The combined output image is sometimes called a “master” of whatever frames you combined, for example, a “master flat” or “master bias”. The task for combining the bias data is **zerocombine** with the altered input parameters and expected output shown in Figure 4. Notice that we will use an average combination operation and identify the object type as “zero” which tells IRAF to collect the bias frames using their header information. Running the task will create a single average combined biased frame titled “Zero.fits” as defined as the output parameter. Reference the bottom frame of Figure 4 to see the combined “master bias”. Creating a master frame will reduce the random noise from the CCD sensor, it will smooth out the pattern of noise providing a better depiction of the bias level that will be subtracted.

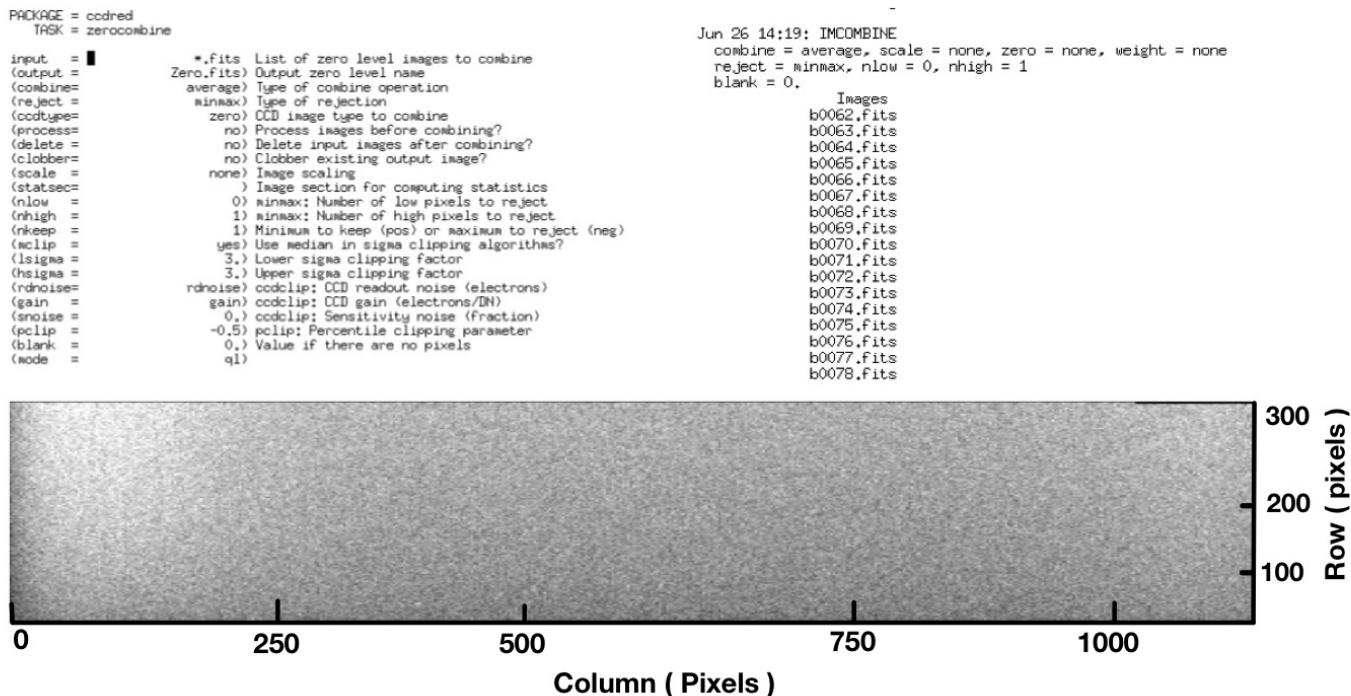

**Figure 4.** Top Left: Initializing the **zerocombine** task with parameters to combine all bias frames by a **minmax** average thus rejecting the minimum and maximum value in each pixel and producing a “master bias” fits file titled Zero.fits. Top Right: The expected output IRAF displays when the combination of selected frames. Bottom: The final combined “master bias” titled “Zero.fits”.

### 3.4. Initializing *ccdproc* with bias correction and trimming the data

The task **ccdproc** is helpful in image reduction and will be run twice. The first run will perform bias subtraction, preflash correction, and trim a defined section of the CCD with the data wanted for analysis. Note that utilizing as much of the CCD is recommended for long-slit data. Therefore, in the **trimsec** parameter, the columns and rows defined as  $[x1:x2,y1:y2]$  will be the portion of the CCD for analysis; any region not defined will be cut off. A benefit of IRAF, is that it provides the option to alter the functions we are trying to fit to our data. When **ccdproc** is run, IRAF will ask in the command line if you want to fit the data interactively. By answering “yes”, you will enter the interactive curve fitting mode where you can change the fitting parameters. Analyzing a low rms and editing the order of the fit is dependent on preference, however, the default settings should be sufficient. Enter “q” to save and exit the fitting mode. If happy with the fit, answer “NO” (in all caps) to the same question in the command line. Once completed, the section defined will be trimmed off and the bias level removed.

```

PACKAGE = ccdred
TASK = Flatcombine

input = ■          *.fits List of flat field images to combine
(output = Flat.fits) Output flat field root name
(combine = average) Type of combine operation
(reject = crreject) Type of rejection
(ccdtype = flat) CCD image type to combine
(process = yes) Process images before combining?
(subsets = yes) Combine images by subset parameter?
(delete = no) Delete input images after combining?
(clobber = no) Clobber existing output image?
(scale = mode) Image scaling
(statsec = ) Image section for computing statistics
(nlow = 1) minmax: Number of low pixels to reject
(nhigh = 1) minmax: Number of high pixels to reject
(nkeep = 1) Minimum to keep (pos) or maximum to reject (neg)
(ncclip = yes) Use median in sigma clipping algorithms?
(lsigma = 3.) Lower sigma clipping factor
(hsigma = 3.) Upper sigma clipping factor
(rdnoise = ccdclp; CCD readout noise (electrons)
(gain = ccdclp; CCD gain (electrons/DN)
(snoise = 0.) ccdclp; Sensitivity noise (fraction)
(pclip = -0.5) pclip; Percentile clipping parameter
(blank = 1) Value if there are no pixels
(mode = q1)

longslit> flatcombine
List of flat field images to combine (*.fits):


Jul 3 13:44: IMCOMBINE
combine = average, scale = mode, zero = none, weight = none
reject = crreject, nclip = yes, nkeep = 1
rdnoise = rdnoise, gain = gain, snoise = 0., hsigma = 3,
blank = 1.

  Images      N  Exp  Mode  Rdnoise  Gain  Scale
b0001.fits   1 900.0 26677. 3. 1.5 0.946
b0002.fits   1 900.0 27315. 3. 1.5 0.924
b0003.fits   1 900.0 26794. 3. 1.5 0.942
b0004.fits   1 900.0 27237. 3. 1.5 0.927
b0005.fits   1 900.0 27124. 3. 1.5 0.931
b0006.fits   1 900.0 27254. 3. 1.5 0.926
b0007.fits   1 900.0 26892. 3. 1.5 0.939
b0008.fits   1 900.0 26922. 3. 1.5 0.938
b0009.fits   1 900.0 26590. 3. 1.5 0.950
b0010.fits   1 900.0 27164. 3. 1.5 0.929
b0011.fits   1 900.0 27208. 3. 1.5 0.928
b0012.fits   1 900.0 26992. 3. 1.5 0.935
b0013.fits   1 900.0 27142. 3. 1.5 0.930
b0014.fits   1 900.0 27399. 3. 1.5 0.921
rFlat.fits  14 900.1 1.1251 3. 1.5 22439.7

, Output image = Flat.fits, ncombine = 28

```

**Figure 6.** Similar to the “Zero.fits” master frame, a “master flat” is created. Left: rejecting any high pixel values caused by the various cosmic rays that fall on the CCD during exposures. Right: Expected output of the combined flat field frames, creating the master frame title “Flat.fits”.

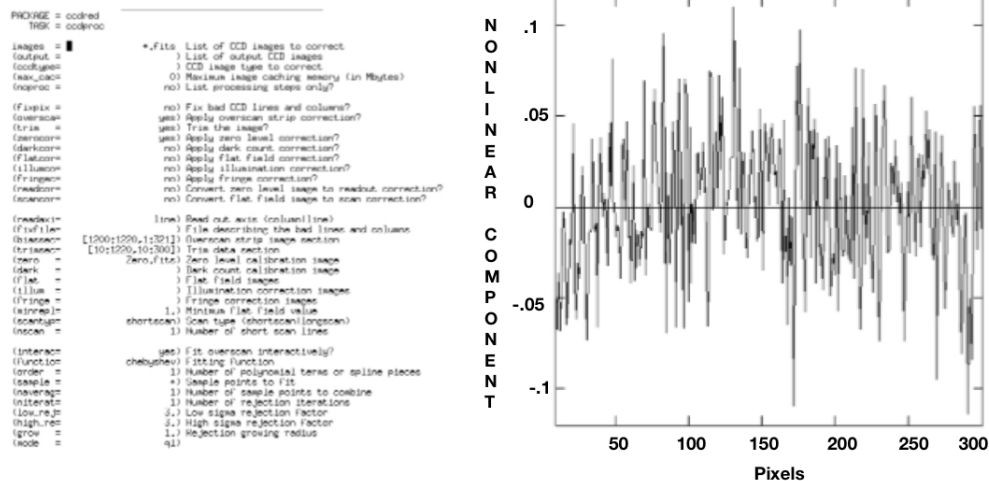


**Figure 5.** Left: Initializing the `ccdproc` with default parameters will create a fit to the bias level of first or second order. Right: An example of the bias level across “Zero.fits”.

### 3.5. Initializing the `fixpix` procedure

Corrupted pixels must be interpolated over with the use of a mask created from the combination of multiple flat field exposures. It is easiest to identify these regions from the flat data. `fixpix` will be the final task that will do interpolate over the bad pixels but needs the combined flat field data. The task for combining flat field data is `flatcombine` with the altered input parameters and expected output shown in Figure 6. Changing the `reject` option to `crreject` will correct for potential radiation events like cosmic rays. This new “master flat” will be used to create the bad pixel mask by entering its output name “Flat” as the input of `ccdmask` task. The default parameters for `ccdmask` do not need to be changed. Once the mask is created, it should be stored in a separate folder so that it does not get altered by the upcoming tasks. The mask will be used as the final step of the reduction process.

### 3.6. Flat correction with `response`

To perform the flat field correction, a “master flat” must be created again. Flat field correction for spectra data has an additional step when normalizing the “master flat” that imaging does not need. Therefore, a function will be fitted to the data along the wavelength axis. To access the task needed to perform the removal and normalization, call the

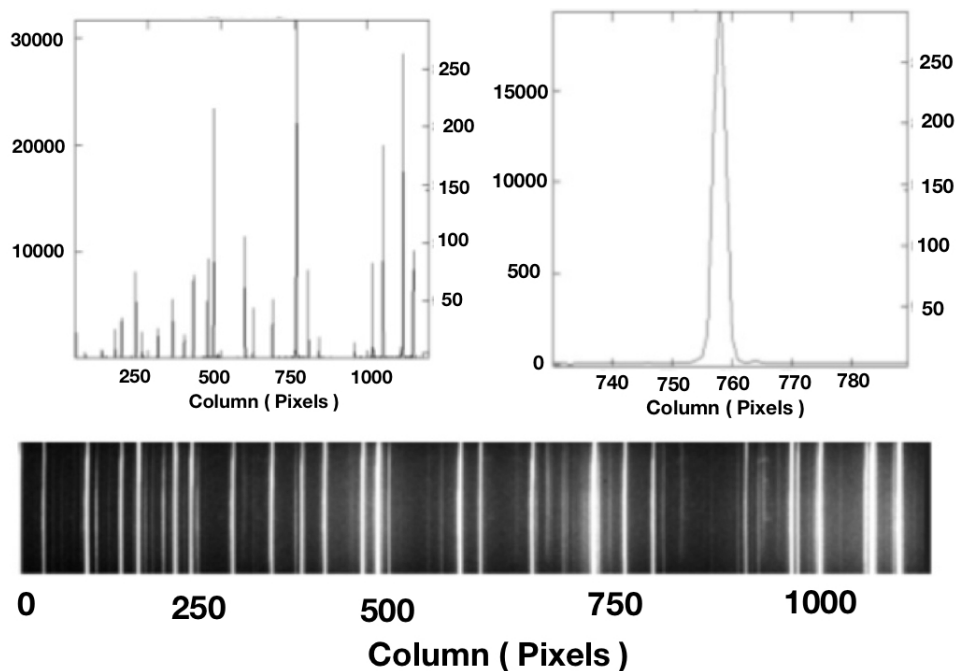

**Figure 7.** Top Left: The emission lineplot that displays the comparison spectra of the HeNeAr lamp used for wavelength calibration. Top Right: Narrowing in on a single line lets the FWHM be defined for the use of an input parameter for future tasks performed for the calibration process. Bottom: The two-dimensional image of the data from the discharge lamp.

primary package **noao** in the command line, then to access the additional packages needed, enter **twodspec**, then subsequently after **long-slit**. The **response** task and default parameters will produce an output image that will be the ratio of the normalized flat frame to the fit created. The function of the fit can be changed to a cubic spline of order 6 by entering **:funct spline3** and **:order 6**. Now the **ccdproc** task can be run once again ( taking into account that it was already run once for the overscan and bias corrections) to perform the flat-field division by entering the output image created with **response** as the **flatcor** input parameter.

### 3.7. Final Correction with *fixpix*

The final correction in the image reduction will be to interpolate over the bad pixels on the CCD. To do so, access the bad pixel mask initially created and stored in a separate folder. Enter it as the input for the **fixpix** task along with the set default parameters. This task will take a few seconds to run, and once completed, a new command line will start, thus concluding the image reduction procedure.


## 4. INTRODUCTION TO WAVELENGTH CALIBRATION

Wavelength calibration is a process that places the reduced two-dimensional spectra on a linear wavelength scale. A fit is created from a comparison spectrum with laboratory wavelengths and is applied to the reduced two-dimensional spectra. The process will be discussed with the IRAF tasks from the **twod.longslit** package.

### 4.1. Arc Frames

An arc or comparison frame is produced by the illumination of a gas-discharge lamp where the light penetrates ionized gas of various elements, thus creating well-characterized emission line spectra. These frames allow for identifying known wavelength measurements as reference for the two-dimensional data. For the long-slit data presented in the left panel of Figure 7, a Helium, Neon, and Argon lamp (HeNeAr) was used at the Bok telescope to produce the comparison spectra. By narrowing in on a single line like on the right of Figure 7, one can help define the line full-width half max (FWHM) for the upcoming wavelength calibration procedure. Analogous to the two-dimensional data shown in Figure 1, the illumination pattern that the HeNeAr lamp produces is shown at the bottom panel.



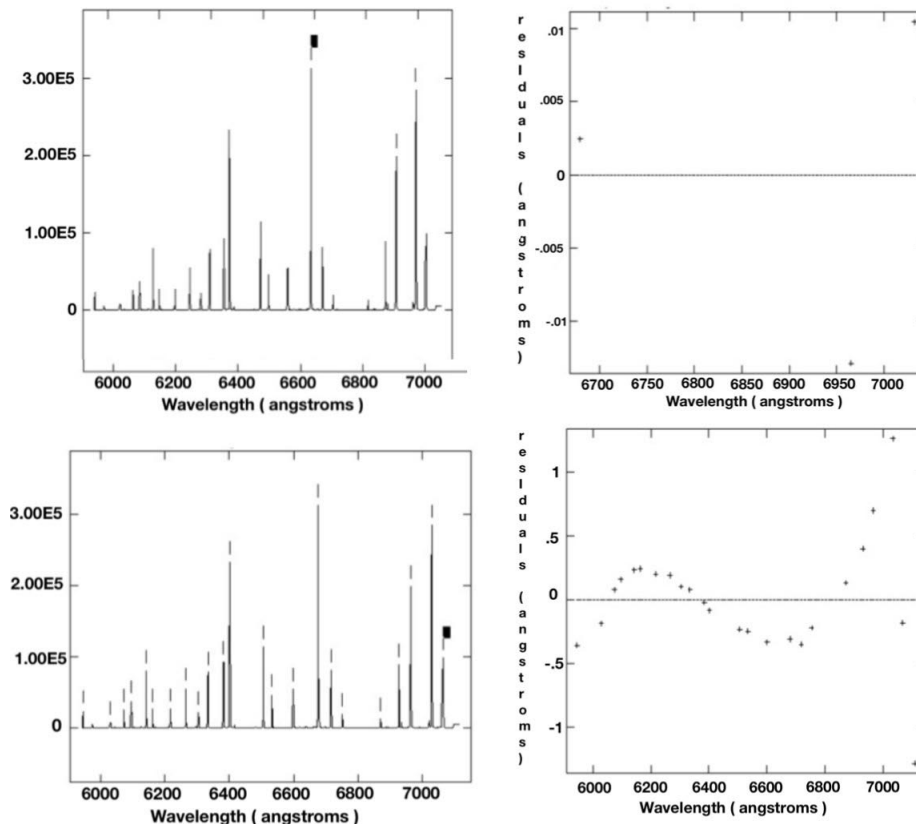

**Figure 9.** Top Left: Initially defining a preliminary fit with a few lines chosen. Top Right: With three lines the preliminary fit is created with a standard ordered fit to let IRAF thus us **identify** determine all of the lines. Bottom right: With all lines defined a better fit is used for the dispersion solution and saved in a database.

be deleted, do so by typing **d** followed by **p** then **f** to have a corrected fit. Yet again, once happy with the overall fit per axis versus residual, save the fit by exiting **fitcoords** with the **q** key.

#### 4.4. Coordinate transformation with *transform*

Lastly, projecting the wavelength solution onto the scientific data will be done with the **transform** task and setting the input parameter to the scientific frame of interest and its corresponding comparison spectra. You may choose whichever name to define as the output frames. Thus, running **transform** will collect the two-dimensional wavelength solution created from the database subdirectory and apply it to the input frames chosen.

#### 4.5. Optional redshift correction with *dopcor*

An additional task can be performed and may be useful for extragalactic spectra. By knowing the systemic redshift of the target, the **dopcor** task can redshift correct the spectra so that any specific emission lines can be shifted to their rest-frame wavelengths.

### 5. CONCLUDING STATEMENTS

An example of the final, reduced 2D spectrum is shown in Figure 12. You may notice that emission skylines are no longer present; a developed Python code removed them. However, the task of background removal and spectra extraction using IRAF can be referenced in *A User's Guide to Reducing Slit Spectra with IRAF* by Phil Massey [2]. Note that skyline removal can be challenging. Thus, it is expected that there will be over-subtraction (darker lines down the location of stronger lines) on the final image. The corrected two-dimensional image can be collapsed into a spatial profile by summing the pixel values along the spatial axis, presented in Figure 12. Figure 13 presents the extracted spectra with a fit (seen in red) on the H $\alpha$  emission line, which is present in the two-dimensional scientific



```

PACKAGE = longslit
TASK = reidentify

reference=      b0057.fits Reference image
images =       b0057.fits Images to be reidentified
(interac=      yes) Interactive fitting?
(section=      middle line) Section to apply to two dimensional images
(newaps =      yes) Reidentify apertures in images not in reference?
(override =    yes) Override previous solutions?
(refit =       yes) Refit coordinate function?

(trace =       yes) Trace reference image?
(step =        10) Step in lines/columns/bands for tracing an image
(nsum =        10) Number of lines/columns/bands to sum
(shift =        0.) Shift to add to reference features (INDEF to sea
(search =       0.) Search radius
(nlost =       0) Maximum number of features which may be lost

(cradius=      10.) Centering radius
(thresho=     10.) Feature threshold for centering
(addfeat=      no) Add features from a line list?
(coordli=     linelistsfidheneat.dat) User coordinate list
(match =      10.) Coordinate list matching limit
(maxfeat=     50) Maximum number of features for automatic identifi
(winssep =     2.) Minimum pixel separation


(database=     database) Database
(logfile=      logfile) List of log files
(plotfile=     ) Plot file for residuals
(verbose=      yes) Verbose output?
(graphic=      stdgraph) Graphics output device
(cursor =     ) Graphics cursor input

answer =       NO Fit dispersion function interactively?
crval =       Approximate coordinate (at reference pixel)
cdelt =       Approximate dispersion
(aidparms=    ) Automatic identification algorithm parameters
(node =       ql)

Reference image = b0057.fits, New image = b0057, Refit = yes
Image Data Found Fit Pix Shift User Shift Z Shift RMS
b0057[*136] 24/24 24/24 -0.105 -0.0976 -1.4E-5 0.243
Fit dispersion function interactively? (no=yes)(NO=YES) (NO): NO
b0057[*136] 24/24 24/24 -0.105 -0.0976 -1.4E-5 0.243
b0057[*126] 24/24 24/24 0.0927 0.0886 1.24E-5 0.222
b0057[*116] 24/24 24/24 -0.0713 -0.0675 -9.8E-6 0.203
b0057[*106] 24/24 24/24 -0.0319 -0.0305 -4.5E-6 0.218
b0057[*96] 24/24 24/24 0.0227 0.0218 3.25E-6 0.207
b0057[*86] 23/24 23/23 0.0447 0.0434 6.66E-6 0.207


```

**Figure 10.** Top: Default parameters of the **reidentify** task that will collect the preliminary fit created to double-check lines. Bottom: the output information indicating the RMS values for each line fits created.

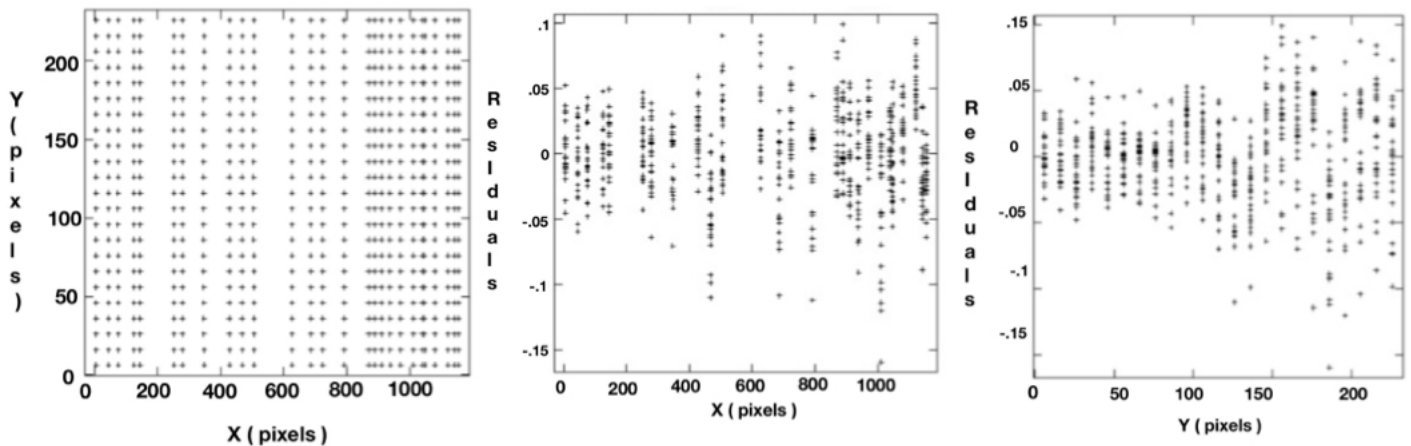


**Figure 11.** Inspecting each fit along the different axes with various key sequences will make sure that the fit has a low RMS for each sequence respectively following left to right along the three panels, (x,y), (x,r), and (y,r).

data. The spectra presented shows the  $H\alpha$  line that's present in an aperture that's placed along the slit of the galaxy disk.

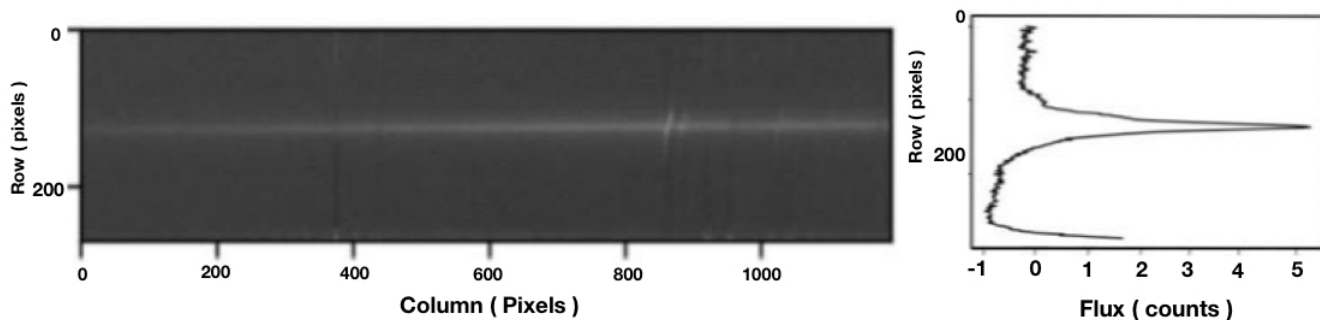

**Figure 12.** What to expect for spectra analysis of a, Left: fully reduced and wavelength corrected two-dimensional spectra of a galaxy with  $H\alpha$  emission lines. Right: An example of analysis regarding the spatial profile of the galaxy.

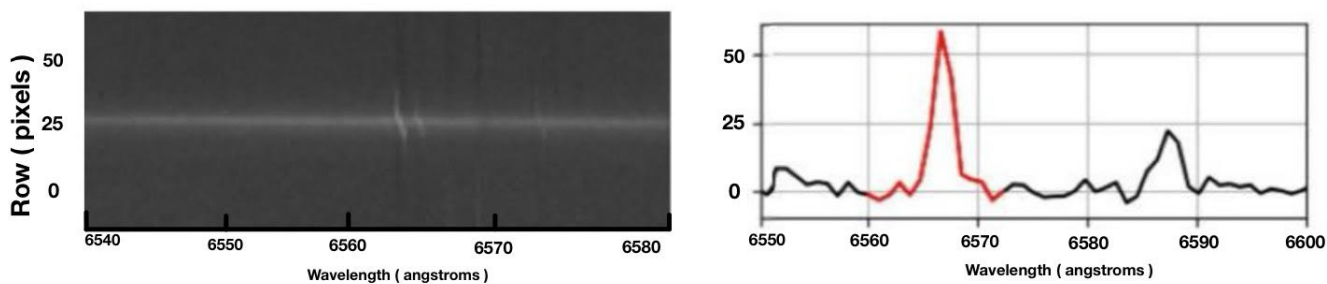

**Figure 13.** The narrowed region on the  $H\alpha$  region of the two-dimensional scientific data and the extracted spectra. In red is a fit specific to the  $H\alpha$  line

## 6. ACKNOWLEDGMENTS

Sara Holeman and Professor Sanchayeeta Borthakur acknowledge the support for this work through the National Science Foundation grants 2108159 and 2408050. Sara Holeman thanks the NSF support through the REU program associated with the grants mentioned-above.

We thank Professor Paul Smith for his input during observations at the Bok telescopes. We also thank the staff at Bok telescope, Steward Observatory for their help and support. We also thank Chris Dupuis, Alejandro Olvera, Mansi Padave, Jackie Monkiewicz, and Brad Koplitz for their advice on observations and comments on this manuscript.

## 7. REFERENCES

- [1] Massey, Phillip. (1997). A User's Guide to CCD Reductions with IRAF. Retrieved 2024, from <https://iraf-community.github.io/doc/iraf.pdf>.
- [2] Massey, P., Valdes, F., & Barnes, J. (1992, April 15). A user's guide to reducing slit spectra with IRAF.