

Sampling of Atmospheric Aerosols by Electrostatic Precipitation for Direct Analyses

Gerd Hermann^{1*}, Georg Lasnitschka¹, Rudolf Matz¹, Alexander Trenin¹ and Walter Moritz²

¹ I. Physikalisches Institut, Justus-Liebig-Universität Gießen, Heinrich-Buff-Ring 16, D-35392 Giessen, Germany

² Fa. Grün-AMS GmbH, D-35630 Ehringshausen, Germany

*Corresponding author (gerd.hermann@exp1.physik.uni-giessen.de)

Abstract. A novel system for aerosol sampling by electrostatic precipitation using graphite platforms as sample collector is presented. Employing standard platforms for commercial analytical instruments, the conception allows fast solid sampling direct element analysis with ETAAS, ETV-ICP-MS/OES, and ETA-CFS without any wet digestive pre-treatment. Other advantages are: highly efficient electrostatic particle collection (>99% for $d = 10 \text{ nm} - 10 \text{ }\mu\text{m}$), reusable sample collectors, omission of filters and chemical reagents. On this basis, an electrostatic precipitator is constructed aiming at a small, relatively uncomplicated instrument. Ten precipitators are arranged in a multi-sampling apparatus for outdoor operation, which simultaneously collect ten samples on same or different collectors for instrumental element analyses, or for microscopic investigations of the collected particles. The precipitator is tested with different model aerosols as well as with atmospheric sampling. Element analysis is carried out with the above mentioned solid sampling methods in external as well as in the authors laboratories.

Keywords: Aerosol analysis; electrostatic precipitation; corona discharge; element analysis; solid sampling; boat technique; graphite platform; scanning electron microscopy; SEM

Introduction

Physical, chemical and biological properties of aerosols, such as total mass and particle number density, spectrum of sizes, variety of shapes, and chemical composition of the particulates are subjects of atmospheric monitoring. Atmospheric particulates range in size over many orders of magnitude from nano-particles of nucleated clusters up to tens of micrometers. The most common method of collecting aerosols is by filtration [1] on fibrous or porous membrane filters. Precipitation by inertial forces is also obtained with cyclones as well as with single-stage or cascade impactors. For chemical analysis the samples are commonly digested in nitric acid, dissolved and subsequently analyzed in liquid form. For direct analysis, atmospheric aerosols have been collected on graphite substrates in cascade impactors [2]. Jaeschke et al. [3] proposed an apparatus for online detection of organic and elemental carbon by non-oxidizing and oxidizing thermal evaporation and subsequent detecting the infrared optical absorption spectrometry of the obtained CO_2 .

However, all approaches of inertial precipitation show cut-off against small aerodynamic diameters. First work with electrostatic aerosol sampling has been carried out by Torge [4] by employing graphite tubes for analysis with atomic absorption spectrometry (AAS) [5]. Measurements by Torge and Hartkamp [6] showed high precipitation efficiency (>98%). Problems of calibration were due to the fact that solid sampling AAS (SSAAS) was not established, at this time.

The first proposal to introduce solid samples directly with the platform, owing to the stationary temperature platform furnace (STPF) concept, as sample boat into the AAS instrument was made by Kurfürst and Rues [7] by employing a SM1 spectrometer (Grün Analytische Meßsysteme GmbH, industry partner of this project). In the interim the method of direct solid sampling has been established [8]. State of the art is the solid sampling boat technique with electrothermal atomization (ETA) AAS (ETAAS) as well as with electrothermal vaporization (ETV) [9,10], but also with Continuum Source Coherent Forward Scattering Spectrometers (CS-CFS) [11,12]. With ETV approaches the sample is converted into an aerosol, which is lead into the analytical instrument, such as ICP or flame. Thereby, the laborious and time-consuming wet-digestive pretreatment can be economized and digestive reagents, that need to be disposed later, are avoided. The graphite boats are reusable and can be cleaned by a simple heating step.

A major problem of sampling atmospheric aerosols onto graphite platforms via inertial impaction [1,2] is the poor efficiency for smaller particles. Even micro-orifice impactors, do not collect particles with cut-off diameters much below $0.1 \text{ }\mu\text{m}$. With previous work [13,14] the authors showed, that particulates in a wide diameter range can be quantitatively collected by electrostatic precipitation (EP) on conductive platforms, such as used in STPF AAS instruments, advantageously those dedicated for the SSAAS (e.g. 'SM 30', Grün AMS GmbH, and 'AAS Solid', Analytik Jena AG) and for ETV for which technique the furnace of SM 30, made by Grün AMS GmbH, is often used [6]. With electrostatic precipitation on such miniature platforms, efficiencies of 99.7% for a particle spectrum generated with laser ablation that

covered a range 5–300 nm were found by measuring with two cascaded precipitators that collected test aerosol particles on the second stage precipitator with and without exciting the corona of the first stage. This is not a measurement of the particle number but of mass collection efficiencies. However, microscopic investigations showed that the test aerosol contained a high fraction of smaller particles. The well-known theory of corona discharges shows that particles with diameters below some 100 nm are mainly charged by diffusion of O_2 -ions, while the larger particles are charged by ions drifted in the electric field [1]. By exposing the aerosol to high ion densities with 10^8 – 10^9 ions s/cm^3 as in our experiments, even nano-particles are charged with high probability. Thus, the electric mobility has a minimum at some 100 nm and increases from this minimum against smaller and against larger diameters as well.

Objective of this work is to develop a sampler for collecting atmospheric aerosols with a broad spectrum of particle diameters on the surfaces of miniature substrates for subsequent analysis which offers the conditions for fast and uncomplicated element analysis by direct solid sampling without any need of chemical digestion and preparation procedures. It should enable to use common analytical instrumentation. The small size of the sample carrier shall allow to construct miniaturized sampling arrangements, well suited for automated instrumentation, and for multi-samplers for simultaneously obtaining sample multitude for analyses with different methods and for comparative measurements.

The approach of electrostatic particle collection shall overcome the problem of poor precipitation efficiencies for smaller particle diameters and of realizing precipitators with low cut-off diameters. As substrates, graphite boats for SSAAS or for ETV-ICP optical emission (ETV-ICP-OES) or mass spectrometry (ETV-ICP-MS), plates cut from silicon semiconductor wafers, highly oriented pyrolytic carbon (HOPC), and plates of glassy carbon (GC) were used for scanning electron microscopy (SEM), glass plates with conductive indium-titanium-oxide (ITO) for optical microscopy.

1 Experimental

The system is presented on Fig. 1. The scheme shows the principle of the EP, and how and to which instruments the sample can be supplied. For ETAAS, ETA-CFS, ETV-ICP the sample boats are introduced directly by manual or automated means, or the samples are vaporized in an ETV that

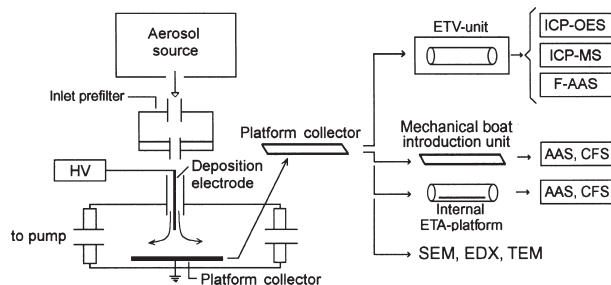


Fig. 1: Scheme of the electrostatic sampling and analyses

first removes the major part of the matrix by pyrolysis, then generates a new aerosol which is fed directly with the argon flow into a flame or the plasma of an ICP. It also can be electrostatically re-precipitated on the internal platform of an ETAAS, or ETA-CFS instrument.

The electrostatic precipitator is designed as point-to-plane configuration. A corona discharge is excited from a sharply grinded tip of a tungsten wire (diameter 120 μm) towards the graphite platform (Fig. 2a). The aerosol flow passes the region of high charge density of the corona, ensuring highly efficient particle charging. A negative voltage of 1.2–1.5 kV is supplied to the tip electrode via a resistor of 30 M Ω causing a mean current of 10–15 μA . The inner diameter of the jet nozzle is 0.5 mm. Its outlet is positioned about 2 mm over and directed against the center of the platform surface. Both the glass capillary and the electrode are adjustable in height. All parts with the electrical connections and O-ring sealings are in an acrylic body (Fig. 2b).

Different commercially available graphite platforms (Fig. 2d) are compatible with matched adapter masks. First the electrical contact to ground potential and current lead-through was realized by thin gold wires, then, with the later improvements, by using a graphite plate shown in Fig. 2a. For microscopic investigations, adaptors are provided for suitable substrates (GC, silicon, HOPC) for SEM.

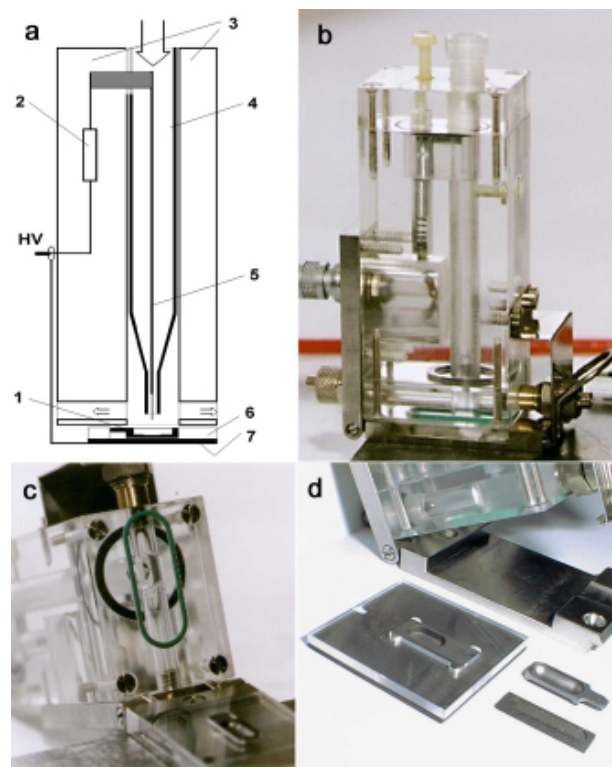


Fig. 2: Electrostatic precipitator. a) Schematic: graphite platform (1), resistor (2), acryl body (3), glass tube with capillary (4), electrode (5), adaptor mask (6), grounded graphite plate (7); photographs: b) complete precipitator; c) opened; d) adaptor mask and commercially available graphite platforms



Fig. 3: a) Ten-fold sampler a) with inlet elutriator and pump; b) front cover opened, ten-fold precipitator with input flow splitter, dosing valves and electric components

The aerosol sampler is operated either by blowing in the model aerosol or by sucking at the outlet ports with a programmable flow-controlled pump (model AirCheck 2000, distributor: MTC GmbH, Mühlheim, Germany). The dosing valves were adjusted with a flow meter (Gilibrator 2, Sensidyne, Clearwater, FL 33760, USA). For multiple sampling, ten of the described precipitators are arranged in a row, connected to a highly symmetric, quasi-isokinetic inlet flow splitter. The flow of each precipitator is adjustable by a dosing valve in each of the outlet flows. The arrangement is mounted in a thermostatically controlled box (Fig. 3) for outdoors measurements. The inlet is through an elutriator (PM5) before entering the flow splitter inside the box. All electrical devices are battery operated, the high voltage was generated by a DC-DC-converter (model 4A12-N4-F, Ultra Volt Inc., Ronkonkoma, NY 11779, USA).

2 Model Aerosols for Laboratory Investigations

For testing of the sampler and its properties, model aerosols had to be established. For aerodynamic particle sizes $d > 0.5 \mu\text{m}$ grinded powders were dispersed with a homemade aerosol gen-

erator. Aerosols with particles sizes $5 \text{ nm} - 1 \mu\text{m}$ were generated by re-condensation of electrothermally vaporized matter.

2.1 Aerosols of larger particle sizes

For testing with larger particle sizes ($100 \text{ nm} - 10 \mu\text{m}$) a laboratory constructed disperser for powders was employed. The powders were grinded with an agate mortar and dispersed with a small magnetically driven stirrer. Subsequently the adequate particle sizes are separated in an elutriator with a low-flow non-turbulent vertical air current. The arrangement is shown in Fig. 4a. Optical microscopy has been carried out with aerosols of sodiumsalicylate precipitated on ITO, scanning electron microscopy (SEM) with precipitations of ash aerosols. For these micrographs substrates of highly oriented pyrolytic carbon (HOPC) have been used. Ash particles with about $100 - 200 \text{ nm}$ in diameter are shown in Fig. 4b).



Fig. 4: a) aerosol generator dispersing particulates; b) ash particles dispersed with generator (a) and electrostatically precipitated on HOPC