

Detection of krypton in xenon for dark matter applications

A. Dobi, C. Davis, C. Hall, T. Langford, S. Slutsky, Y.-R. Yen

Department of Physics, University of Maryland, College Park MD, 20742 USA

Abstract

We extend our technique for observing very small concentrations of impurities in xenon gas to the problem of krypton detection. We use a conventional mass spectrometer to identify the krypton content of the xenon, but we improve the sensitivity of the device by about five orders of magnitude with a liquid nitrogen cold trap. We find that the absolute krypton concentration of the xenon can be inferred from the mass spectrometry measurements, and we identify krypton signals at concentrations as low as 0.5×10^{-12} mol/mol (Kr/Xe). This technique simplifies the monitoring of krypton backgrounds for WIMP dark matter searches in liquid xenon.

Key words:

xenon, krypton, cold trap, mass spectrometry, dark matter

1. Introduction

One of the primary challenges faced by liquid xenon WIMP dark matter experiments is the presence of trace amounts of radioactive krypton. Xenon itself has no long-lived radioactive isotopes which might act as background sources, but krypton includes the troublesome anthropogenic isotope ^{85}Kr , a beta emitter with a Q value of 687 keV and a half-life of 10.76 years. ^{85}Kr is created in nuclear power plants and released into the earth's atmosphere during fuel reprocessing, and its isotopic fraction at present is about 2×10^{-11} mol/mol ($^{85}\text{Kr}/^{nat}\text{Kr}$) [1]. Xenon, on the other hand, is extracted from the atmosphere with a residual krypton concentration typically ranging from 10^{-9} to 10^{-6} mol/mol ($^{nat}\text{Kr}/\text{Xe}$). Although this implies that the absolute concentration of ^{85}Kr in xenon is rather small, the ^{85}Kr beta decay is nevertheless highly problematic for dark matter experiments because these

decays are not suppressed by self-shielding and because krypton cannot be separated from xenon with conventional chemical purifiers.

Modern dark matter experiments such as XMASS[2], XENON100[3, 4], and LUX[5] are expected to probe WIMP cross-sections down to $\sim 7 \times 10^{-46}$ cm² for a 100 GeV WIMP mass. This requires that the residual krypton concentration of the xenon target material be no more than $\sim 3 \times 10^{-12}$ mol/mol ($^{nat}\text{Kr}/\text{Xe}$)¹ so that ^{85}Kr is a sub-dominant background source. Future xenon experiments intend to probe WIMP cross sections down to 10^{-48} cm², which will require reducing the krypton concentration by another factor of sixty[5, 6, 7].

To achieve these ultra-low krypton concentrations, commercially procured xenon must undergo additional processing via distillation[8] or gas chromatography[9]. Once this processing is complete, the residual krypton content of the xenon can be analyzed via low background counting [8], mass spectrometry aided by chromatography [9], or atomic trap trace analysis (ATTA) [10]. The low background counting and chromatographic methods have achieved sensitivities of about $\sim 1 \times 10^{-12}$ mol/mol, while the proposed ATTA technique is expected to achieve a sensitivity of 3×10^{-14} mol/mol. These monitoring techniques can confirm that the krypton removal processing has been successful prior to full detector operations, and can constrain the total background count rate due to krypton in the WIMP search data.

In this article we show that very small krypton concentrations can be observed in xenon gas using a mass spectrometry technique which we previously developed to detect electronegative impurities in xenon[11]. The method is simple, highly sensitive, and easily mastered, and it could be quickly adopted and applied by many working dark matter experiments.

2. Xenon cold trap mass spectrometry

We use a residual gas analyzer (RGA) mass spectrometer to analyze our xenon by introducing a small quantity of the gas into the RGA's vacuum enclosure through a leak valve. Since the partial pressure of each component species is proportional to both its absolute concentration and to the flow rate through the analysis system, by controlling for the flow rate the partial pressures can be interpreted in terms of the absolute concentrations. The

¹All concentrations in this article refer to the natural krypton to xenon ratio, measured in units of mol/mol, unless otherwise indicated.

measurement is calibrated by preparing samples of xenon gas with known impurity concentrations for the various species of interest by directly mixing known quantities of impurities with a known amount of xenon. The flow rate can be controlled either by measuring the actual flow rate in real time or by using a standard leak-valve setting whose flow rate was previously calibrated.

Since the partial pressures of all species are proportional to the flow rate, the RGA signals can be vastly increased simply by opening the leak valve further. However, the RGA cannot be operated above some maximum total pressure, typically about 10^{-5} Torr, and the total pressure is dominated by the xenon present in the gas sample. This limits the maximum flow rate that can be used. For example, if krypton can be detected by the RGA at a partial pressure of $\sim 10^{-12}$ Torr, and the xenon pressure is 10^{-5} Torr, then the limit of detection is about one part in 10^7 . Since we are interested in krypton concentrations at the level of 10^{-12} mol/mol, this is inadequate for our purposes.

We solve the saturation problem simply by removing most of the xenon from the gas sample with a liquid nitrogen cold trap placed between the leak valve and the RGA. The cold trap holds the xenon pressure constant at the vapor pressure of xenon ice, which allows the flow rate to be vastly increased without saturating the RGA. In our previous paper, we showed that impurity species such as oxygen, nitrogen, and methane pass through the cold trap in large quantities, and that their partial pressures, corrected for flow rate, remain proportional to their absolute concentrations. The sensitivity to oxygen, nitrogen, and methane was found to be 0.66×10^{-9} , 9.4×10^{-9} , and 0.49×10^{-9} (mol/mol), respectively [11].

Here we extend the technique to observe krypton in xenon. We expect that krypton could be observed in very small quantities by the RGA because there are very few background species which could obscure the krypton signal. In fact, we find that we are able to detect krypton in xenon at a concentration of 0.5×10^{-12} mol/mol, which makes this technique better than or comparable to existing methods, and sensitive enough to be useful for working dark matter experiments.

3. Apparatus and procedures

A diagram of our analysis apparatus is shown in Figure 1. The xenon of interest is admitted into the analysis system through an ultra-high vacuum leak valve (Kurt Lesker part number VZLVM940R). It passes through a

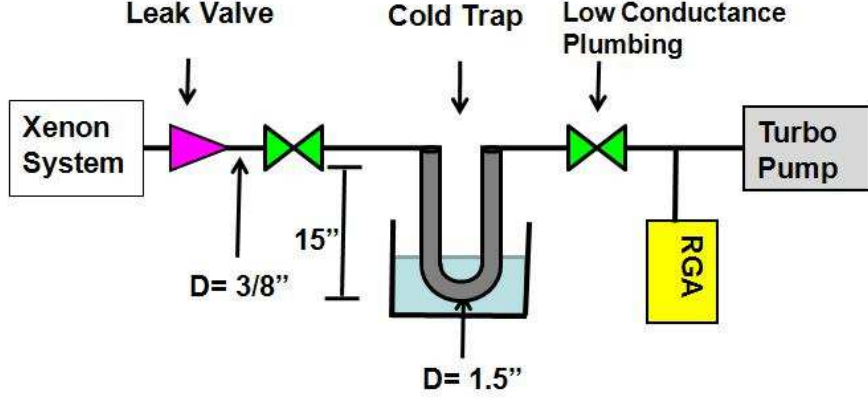


Figure 1: Diagram of the cold trap mass spectrometer analysis system.

liquid nitrogen cold trap and a section of low-conductance plumbing before reaching an SRS RGA200 mass spectrometer. The low-conductance plumbing is necessary to reduce the xenon partial pressure from 1.8×10^{-3} Torr (its vapor pressure at liquid nitrogen temperature [12]) to $< 10^{-5}$ Torr. This insures that the RGA remains unsaturated. We use a fully open hand valve for the low conductance element.

The cold trap is constructed from 1.5" OD stainless-steel tubing, welded into a U shape with a radius of 2.5" inches and 12.5" linear inlet and outlet legs². The tubing diameter is chosen to allow a significant amount of xenon to be analyzed before the growth of xenon ice blocks the flow of gas through the analysis system. We find that the liquid nitrogen level must be high enough to submerge the bottom of the cold trap U, but otherwise its level is not critical.

We use several procedures to calibrate and monitor the flow rate through the leak valve. First, we calibrate the flow rate directly for a variety of leak valve settings using a fixed volume of xenon gas (~ 1 liter) at the leak valve input. We measure the pressure drop in this volume with a capacitive manometer as the gas flows through the leak valve to infer the flow rate for each leak valve setting. We find that the flow rate is repeatable to within 10% simply by returning the leak valve to the same indicator marking on its dial.

²We have also constructed working cold traps from standard vacuum plumbing components with 2.75" CF flanges.

Secondly, we also measure the flow rate directly using a MKS model 179A mass flow meter which has been calibrated for use with xenon gas. Third, when the xenon under analysis contains a small, constant concentration of a tracer gas which is unaffected by the cold trap, such as argon or helium, then the partial pressure of the tracer can be used to accurately monitor the leak rate in real time with the RGA itself. This eliminates the systematic error due to the leak valve dial setting and RGA gain drift. We use this method in Section 4.

To calibrate the partial pressure measurements of the RGA in terms of the true krypton concentration, we insert known quantities of krypton into the xenon using a krypton gas cylinder (99.999% krypton purity). The injection volume is 13.8 ccs of plumbing monitored by a pressure gauge and isolated by two valves. This section of plumbing is pumped out to vacuum, and then re-pressurized with krypton from the krypton cylinder. The krypton is combined with the xenon by flowing the xenon through the injection volume and collecting the gases in a recovery bottle where they mix.

To perform a measurement, first we submerge the cold trap in liquid nitrogen while it is pumped to ultra-high vacuum by the turbo-molecular pump. We then open the leak valve in two steps. In the first step, we use a very small leak rate, less than 10^{-4} standard liters per minute (SLPM), which allows xenon ice to form in the cold trap, establishes the fixed xenon partial pressure, and flushes some trace background gases out of the analysis system plumbing. We wait for several minutes for the partial pressures of all species to stabilize, and then we open the leak valve to the desired flow rate for purity analysis. In general, the best sensitivity is obtained by using the maximum possible flow rate. In some cases the flow rate is limited by partial pressure of non-xenon impurity species such as oxygen, nitrogen, or argon. Since these species are not removed by the cold trap, their presence in the xenon gas will eventually cause the RGA to saturate as the flow rate is increased. For the very best sensitivity, the xenon should be free from extraneous impurity species.

4. Response of the analysis system to krypton

In our first series of experiments, we confirm that the krypton partial pressure observed by the RGA is indeed linear in the true concentration by injecting known amounts of krypton into our 2.8 kg xenon supply. For these measurements, the xenon supply bottle continuously feeds xenon into

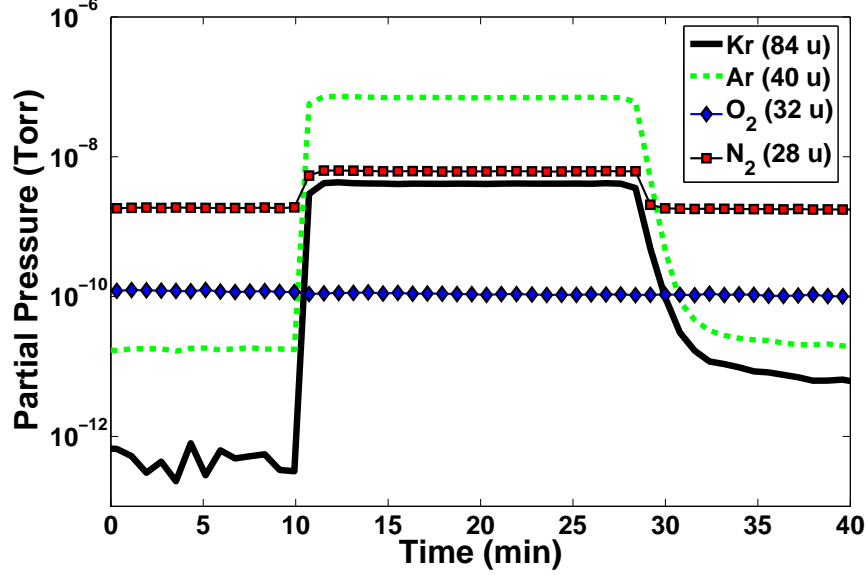


Figure 2: The results for a typical purity measurement with constant flow rates. At $t = 10$ minutes the leak valve is opened and a measurement is made with a flow rate of 0.1 SLPM. This sample of xenon contained 20×10^{-9} mol/mol nitrogen, 10^{-6} mol/mol argon, and 66.4×10^{-9} mol/mol krypton. At $t = 28$ minutes the leak valve is closed.

the system through a regulator, maintaining a constant pressure at the leak valve input. This insures that the leak rate into the analysis system is nearly constant throughout the measurement, which simplifies the data analysis.

To precisely monitor the leak valve flow rate in real time, we use argon as a tracer gas. Our 2.8 kg xenon supply contains an argon concentration of about 10^{-6} mol/mol, and since the argon level is constant from one injection experiment to the next (because the injected gas is 99.999% krypton, with only trace quantities of argon), the argon partial pressure serves as a convenient proxy for the gas flow rate through the analysis system. Under these conditions, the leak valve can be opened to an arbitrary setting, and the krypton-to-argon partial pressure ratio should be proportional to the true krypton concentration.

A typical measurement is shown in Figure 2. At $t = 0$, xenon ice has already been established in the cold trap, and the measurement begins at $t = 10$ minutes. Krypton, argon, and nitrogen are clearly present in the sampled gas, while the oxygen concentration is less than 0.7×10^{-9} mol/mol.

$\Delta\rho(\text{Kr})$ (10^{-9} mol/mol)	P_{Kr} (84 u) (10^{-9} Torr)	P_{Ar} (40 u) (10^{-9} Torr)	P_{Kr}/P_{Ar} (Torr/Torr)
0	4.11	70.3	0.0584
7.37	5.90	91.2	0.0647
	6.02	91.9	0.0655
18.41	7.00	93.7	0.0747
	6.60	87.8	0.0752
33.14	7.88	89.6	0.0879
	7.51	86.0	0.0873
	8.05	92.0	0.0876
	7.76	89.1	0.0871

Table 1: Krypton and argon partial pressures as a function of the injected krypton concentration ($\Delta\rho(\text{Kr})$). Each gas sample was measured at least twice to gauge the systematic uncertainty in the partial pressure measurements for a fixed concentration of krypton. We infer an uncertainty of 5% in the krypton-to-argon partial pressure ratio.

At $t = 28$ minutes the leak valve is closed.

In Figure 3 and Table 1 we show the krypton-to-argon partial pressure ratio as a function of the amount of krypton which we inject into our xenon supply. The ratio is linear in the injected concentration, which confirms that the cold trap allows the krypton to pass through as desired. In total we injected 33.1×10^{-9} mol/mol of krypton, which resulted in a total increase in the krypton-to-argon ratio of a factor of 1.50 relative to the vendor-supplied xenon. From this we infer that the krypton concentration was $(66.4 \pm 4) \times 10^{-9}$ mol/mol before our injections, and $(99.5 \pm 4) \times 10^{-9}$ mol/mol after injections.

5. Detection of krypton at the 10^{-12} mol/mol level

The analysis system response to krypton is 0.79 Torr/(SLPM \cdot mol/mol). Since the fluctuations in the RGA reading at 84 AMU are $\sim 3 \times 10^{-13}$ Torr, we expect that a concentration of 1×10^{-12} mol/mol krypton could be detectable at a flow rate of ~ 0.4 SLPM. However, as shown in Figure 3, our xenon supply contains significant argon and nitrogen. Since these trace gases are not removed by the cold trap, they will cause the RGA to saturate at flow rates above 0.1 SLPM, with the argon being the leading problem. This would limit our krypton sensitivity to about 4×10^{-12} mol/mol. To detect krypton at lower concentrations it is necessary to remove these trace impurities.

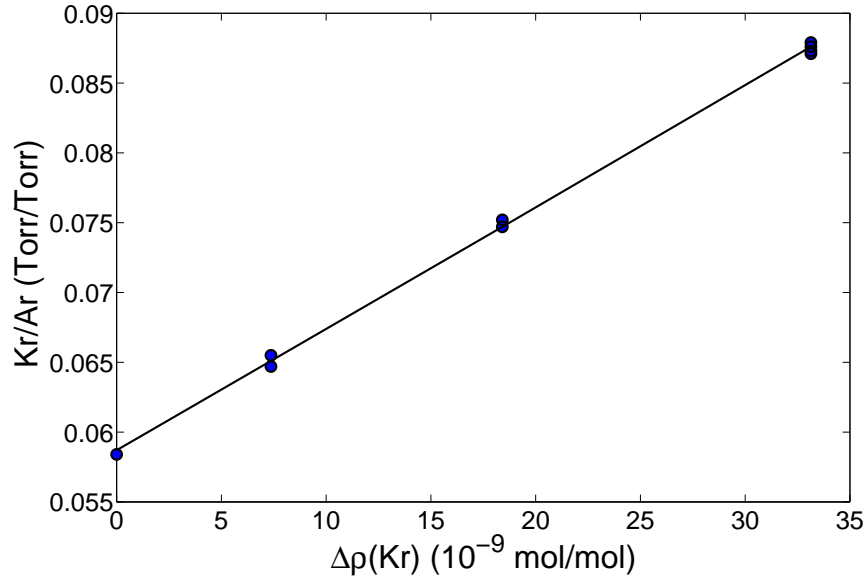


Figure 3: Krypton-to-argon partial pressure ratio versus the injected krypton concentration ($\Delta\rho(\text{Kr})$). Repeated measurements are shown as separate data points. The non-zero y-intercept value is due to the krypton present in our vendor supplied xenon before our injections.

Nitrogen can be removed from xenon using standard getters [13], but argon cannot. We first tried to remove the argon by freezing the xenon in its supply bottle with liquid nitrogen and pumping on the vapor with the turbo-molecular pump. This strategy proved to be inefficient, probably because the argon is trapped in the xenon ice. However, we successfully purified a small quantity of xenon using the cold trap itself, since the cold trap captures xenon while allowing argon and krypton to pass through. We allowed approximately seven standard liters of xenon to slowly leak into the cold trap over a period of two hours and then recovered the xenon into a storage bottle. We repeated this process a total of seven times. The average argon and krypton reduction factors for each pass were 180 and 7.4. This produced a 40 gram sample of xenon with much less than 1×10^{-12} mol/mol argon and krypton.

Starting with this 40 gram sample of de-argonated and de-kryptonated xenon gas, we created xenon with known krypton concentration by mixing it with small quantities of the 2.8 kg xenon supply, which contained $99.5 \times$

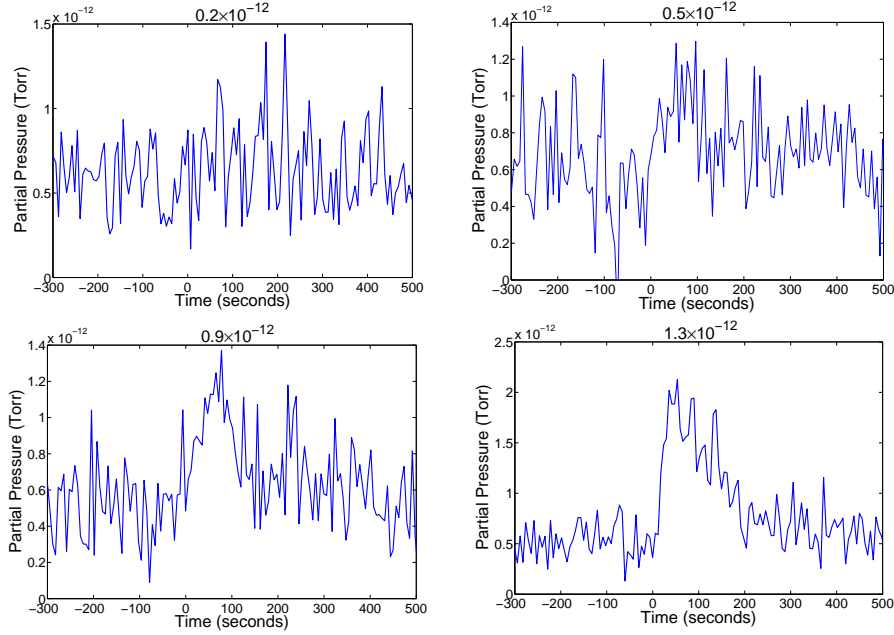


Figure 4: RGA response to the smallest concentrations of krypton. The krypton signals decay in time because the flow rate is decreasing due to the small amount of xenon available for these measurements. Upper left: 0.2×10^{-12} mol/mol. Upper right: 0.5×10^{-12} mol/mol. Lower left: 0.9×10^{-12} mol/mol. Lower right: 1.3×10^{-12} mol/mol.

10^{-9} mol/mol of krypton after the experiments described in Section 4. For example, to achieve 0.5×10^{-12} mol/mol of krypton, we added 0.2 milligrams of our krypton-rich xenon supply to the 40 grams of de-kryptonated xenon. The xenon was then analyzed by the cold trap mass spectrometry technique. Between each run we removed the extra krypton, and a new sample was mixed starting again with de-kryptonated xenon.

RGA partial pressure plots are shown in Figure 4 for xenon with 0.2, 0.5, 0.9, and 1.3×10^{-12} mol/mol of krypton. For the data shown in these plots we open the leak valve to its maximum setting at $t = 0$ seconds. Since the total amount of xenon available for these measurements is modest, the flow rate immediately peaks at 1.5 SLPM and then decreases during the measurement due to the decreasing pressure at the leak valve input. The resulting partial pressure plots follow this pattern, as shown in Figure 4.

Clear krypton signals are seen for concentrations of 0.9×10^{-12} mol/mol and larger, while the 0.5×10^{-12} mol/mol sample gives a marginal signal.

$\rho(\text{Kr})$ (10^{-12} mol/mol)	Avg. P_{Kr} (10^{-12} Torr)	Avg. Flow (SLPM)	Kr Fig. of Merit (10^{-12} Torr/SLPM)	Dev. from Fit (%)
0.5 ± 0.2	0.406	1.32	0.308 ± 0.0154	3.8
0.9 ± 0.2	0.646	1.27	0.507 ± 0.0254	5.1
1.3 ± 0.2	1.15	1.28	0.901 ± 0.0450	-12.4
1.7 ± 0.2	1.37	1.28	1.07 ± 0.0528	1.2
5.1 ± 0.3	3.81	1.30	2.93 ± 0.147	8.6
17.1 ± 0.9	15.7	1.29	12.2 ± 0.608	-13.7
171.1 ± 8.8	127	1.27	99.5 ± 4.97	6.9
* 1711 ± 88	1450	1.31	1110 ± 55.3	-3.5
1711 ± 88	1330	1.29	1030 ± 51.3	4.1

Table 2: Results of krypton detection experiments with 40 grams of highly purified xenon. The various krypton concentrations were created by mixing with xenon containing 99.5×10^{-9} mol/mol krypton, except the sample labeled (*), which was created by injecting 99.999% krypton from a krypton gas cylinder. For prepared samples of 1.7×10^{-12} mol/mol or less the uncertainty in the concentration is dominated by the minimum sensitivity to krypton in the highly purified xenon, which we take to be 0.2×10^{-12} mol/mol. For concentrations above 1.7×10^{-12} the uncertainty in the concentration is 5%, dominated by the uncertainty of the concentration of the krypton-rich xenon supply. The krypton figure of merit is the average partial pressure divided by the average flow rate. The last column shows the deviation from the linear fit shown in Figure 5.

No significant deviation from background is seen for the 0.2×10^{-12} mol/mol sample, which agrees with our expected sensitivity of 0.25×10^{-12} mol/mol (inferred by assuming partial pressure fluctuations of 0.3×10^{-13} Torr and a 1.5 SLPM flow rate.)

To quantify the krypton concentration, we calculate the average partial pressure in a 60 second window around its maximal value after subtracting the background level, and we divide by the average flow rate measured by the MKS mass flow meter. Since the same leak valve setting (fully open) was used in each dataset, the average flow rate varies by less than 2% in all datasets. As shown in Figure 5 and Table 2, the ratio of average pressure to average flow rate is proportional to the true krypton concentration over four orders of magnitude. The largest deviation from the fitted line is 13.7%, and the dataset at 0.5×10^{-12} mol/mol deviates from the fitted line by only 4%. This indicates that the small krypton signal at 0.5×10^{-12} mol/mol is likely to be genuine, and demonstrates sensitivity to krypton at concentrations less than 1×10^{-12} .

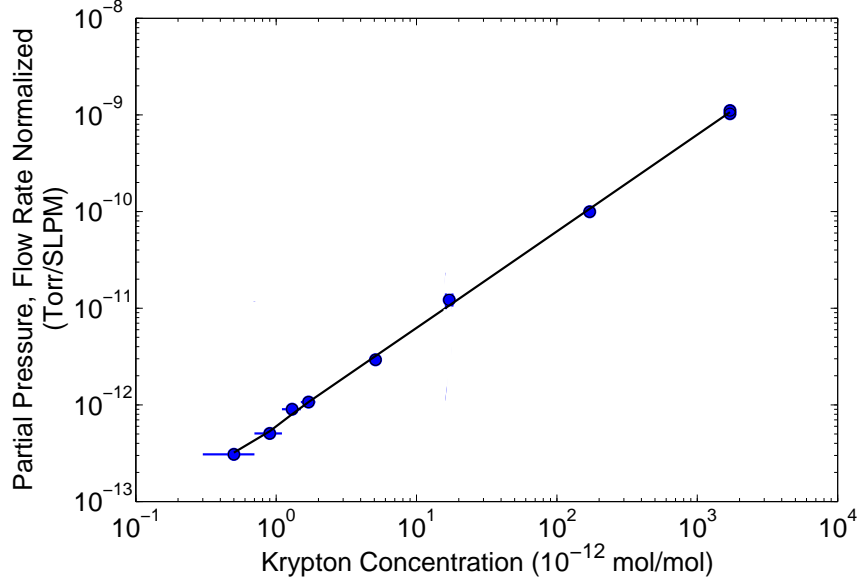


Figure 5: Krypton detection figure-of-merit as a function of true krypton concentration. The krypton partial pressure data (at 84 AMU) is averaged in a 60 second window around its peak value and then normalized to the average flow rate. The maximum deviation from a linear fit over four orders of magnitude is 13.7%

To confirm that the absolute krypton concentration of our experiments is not in error, we performed one final krypton injection from the krypton cylinder into our 40 grams of de-kryptonated xenon at a concentration of 1.7×10^{-9} mol/mol. This dataset gives a krypton figure of merit which agrees with the equivalent dataset produced by mixing to within 8%.

6. Conclusion

We have extended the xenon cold trap mass spectrometry technique to detect trace quantities of krypton in xenon gas. We find that krypton passes through the cold trap largely undisturbed, and that the resulting partial pressure is proportional to the true concentration after accounting for the flow rate. Using this method we have detected krypton concentrations as low as $(0.5 \pm 0.2) \times 10^{-12}$ mol/mol $^{nat}\text{Kr}/\text{Xe}$.

We believe our sensitivity could be significantly improved by using faster flow rates, which could be achieved by using a larger sample of highly pu-

rified xenon. In principle we see no reason why an additional factor of ten improvement could not be achieved by using ten times the amount of xenon, which would make the technique useful for future WIMP dark matter experiments. In any case, the sensitivity demonstrated here will already be useful for krypton monitoring programs at existing detectors.

7. Acknowledgments

This work was supported by the National Science Foundation under award number PHY0810495.

References

- [1] P. Collon, W. Kutschera, Z.-T. Lu, *Ann. Rev. Nucl. Part. Sci.* 54 (2004) 39–67.
- [2] K. Abe, *J. Phys. Conf. Ser.* 120 (2008) 042022.
- [3] E. Aprile, et al., *Phys.Rev.Lett.* 105 (2010) 131302.
- [4] E. Aprile, et al. [arXiv:1101.3866](#).
- [5] D. McKinsey, et al., *J. Phys. Conf. Ser.* 203 (2010) 012026.
- [6] K. Arisaka, et al., *Astroparticle Physics* 31 (2) (2009) 63 – 74.
- [7] L. Baudis [arXiv:1012.4764](#).
- [8] K. Abe, et al., *Astropart. Phys.* 31 (4) (2010) 290–296.
- [9] A. Bolozdynya, et al., *Nucl. Instrum. Meth A* 579 (1) (2007) 50–53.
- [10] D. N. McKinsey, C. Orzel, *Nucl. Instrum. Meth. A* 545 (2005) 524–531.
- [11] D. Leonard, et al., *Nucl. Instr. Meth. A* 621 (1-3) (2010) 678 – 684.
- [12] A. Grutter, J. Shorrocks, *Nature* 204 (1964) 1084–1085.
- [13] A. Dobi, et al., *Nucl. Instr. Meth. A* 620 (2-3) (2010) 594 – 598.