

# Stimulated Raman Spectroscopy with Tunable Visible Broadband Probe Pulse Generated by Kerr-Instability Amplification

Nathan G. Drouillard, TJ Hammond

## Abstract

Femtosecond, broadband stimulated Raman spectroscopy is a popular approach to measuring molecular dynamics with excellent signal-to-noise and spectral resolution. We present a new method for broadband stimulated Raman spectroscopy that employs Kerr-instability amplification to amplify the supercontinuum spectrum from sapphire and create a highly tunable Raman probe spectrum spanning from 530 to 1000 nm ( $-6000$  to  $2800\text{ cm}^{-1}$ ). Our method, called Kerr-instability amplification for broadband stimulated Raman spectroscopy (KAB-SRS) provides an alternative to optical parametric amplifiers by producing a broader and more tunable spectrum at a fraction of the cost. We demonstrate the effectiveness of our method by measuring the stimulated Raman loss spectrum of 1-decanol.

## Keywords

Tunable, broadband femtosecond stimulated Raman spectroscopy, stimulated inverse Raman spectroscopy, anti-Stokes, ultrafast Raman loss spectroscopy, 1-decanol, Kerr-instability amplification.

## Introduction

With molecular vibrations occurring on the timescale of femtoseconds ( $1\text{ fs} = 1 \times 10^{-15}\text{s}$ ) and rotations for small molecules on the order of picoseconds ( $1\text{ ps} = 1 \times 10^{-12}\text{s}$ ), spectroscopic measurements that probe the structure of matter benefit greatly from ultrafast light sources<sup>1</sup>. While there exist ultrafast electronic spectroscopies, these techniques often make it difficult to measure structural dynamics<sup>1,2</sup>. Therefore, vibrational spectroscopy is a more suitable approach to elucidating the desired molecular information.

Most molecules possess several Raman active modes that provide information on molecular vibrations, making Raman spectroscopy a straightforward alternative to infrared (IR) spectroscopy, which requires far-IR light sources and detectors. Especially with regards to aqueous samples that are dominated by infrared absorption, Raman spectroscopy is a viable alternative<sup>3,4</sup>.

Femtosecond stimulated Raman spectroscopy (FSRS) has become a common approach to studying molecules and reaction pathways on ultrafast timescales<sup>5-7</sup>. FSRS offers several advantages over spontaneous Raman, namely its fast acquisition time, superior signal-to-noise

ratio, and fluorescence rejection making it a background-free technique<sup>2,8</sup>. The theory of FSRS has been published in detail elsewhere<sup>1,2</sup>.

Typically performed as a three-pulse experiment, FSRS can be thought of as stimulated Raman spectroscopy (SRS) with an additional pulse called the actinic pump pulse that provides the temporal resolution of the measurement. Without the actinic pump, but while still having a femtosecond, broadband Raman probe pulse, the name broadband stimulated Raman spectroscopy (BSRS) is often used<sup>9,10</sup>. In this work, we are not concerned about the temporal resolution provided by an actinic pulse, but we sought to develop a novel technique for generating the broadband Raman probe pulse, which ultimately defines the bandwidth of detection in BSRS. Our efforts are aligned with recent developments in the field of FSRS which have been focussed on the tunability of ultrafast laser pulses<sup>11, 12</sup>.

A common approach to generating a broadband probe spectrum is supercontinuum generation in sapphire<sup>2</sup>. Downsides to supercontinuum generation include intense residual pump light and less intense light at longer wavenumbers. Alternatively, one may use an optical parametric amplifier (OPA) or its noncollinear counterpart (NOPA)<sup>2,13,14</sup>. The bandwidth of this method is typically limited to  $1200 \text{ cm}^{-1}$  due to gain narrowing<sup>2,15</sup> and is limited to gain media that possess  $\chi^{(2)}$  optical nonlinearities, far less common than those that exhibit a  $\chi^{(3)}$  response<sup>16</sup>. Fang et al have developed an innovative alternative that uses cascaded four-wave mixing (CFWM) to produce multiple sidebands, each having a FWHM bandwidth of up to  $1600 \text{ cm}^{-1}$ , with a total spectral region of detection spanning  $100\text{-}4000 \text{ cm}^{-1}$ <sup>17,18</sup>. Yet still, one must scan using multiple Raman probe beams to span the entire Raman spectrum from  $0\text{-}4000 \text{ cm}^{-1}$ . Conversely, we present a new method that generates a tunable, broadband probe that has a continuous spectrum spanning from  $-6000 \text{ cm}^{-1}$  to  $0 \text{ cm}^{-1}$  on the anti-Stokes side, making it possible to make such measurements with one individual probe spectrum. Furthermore, we can generate a probe spectrum that spans from  $0$  to  $2800 \text{ cm}^{-1}$  on the Stokes side, by changing the filter used prior to amplifying the supercontinuum.

Since we operate with a blue-shifted probe in this work, we detect stimulated inverse Raman scattering (IRS)<sup>19-21</sup>, sometimes referred to as Raman loss. Hence, techniques such as ours are often referred to as stimulated inverse Raman spectroscopy, anti-Stokes (femtosecond) stimulated Raman spectroscopy, or ultrafast Raman loss spectroscopy (URLS)<sup>22-25</sup>. For this work, we will use the convention that IRS refers to the physical process that is measured via URLS<sup>25</sup>. Stimulated inverse Raman scattering is an intensity dependent process where the intensity of the anti-Stokes field grows exponentially with the pump intensity,  $I_{as}(z) \propto \exp(|g|I_{pu}z)$ , with  $g = \frac{2\pi\omega_{pu}}{nc} \text{Im}[\chi_{NL}]$  and where  $\chi_{NL}$  is the nonlinear susceptibility<sup>26,27</sup>.

We present Kerr-instability amplification for broadband stimulated Raman spectroscopy (KAB-SRS), which uses Kerr-instability amplification (KIA) rather than CFWM. KIA extends four-wave mixing to the high-intensity regime, leading to a relaxed phase-matching condition and enabling broadband amplification across nearly an octave of bandwidth. The amplified spectrum from KIA is tunable to both the blue (anti-Stokes) and red (Stokes) sides of the pump wavelength. The theory for KIA and the realization of broadband amplification are outlined in our previous work<sup>28-31</sup>. As established in our prior research, KIA provides a cost-effective, versatile route for generating broadband Raman probe spectra. We can amplify down to  $500 \text{ nm}$  with a  $785 \text{ nm}$  pump ( $-7000 \text{ cm}^{-1}$ ) in one single spectrum, but we are limited by dispersion of

the reimaging optics that come after amplification rather than phase matching in the case of an OPA or NOPA<sup>25</sup>. We show the efficacy of our method by measuring the Raman loss spectrum of 1-decanol.

## Methodology

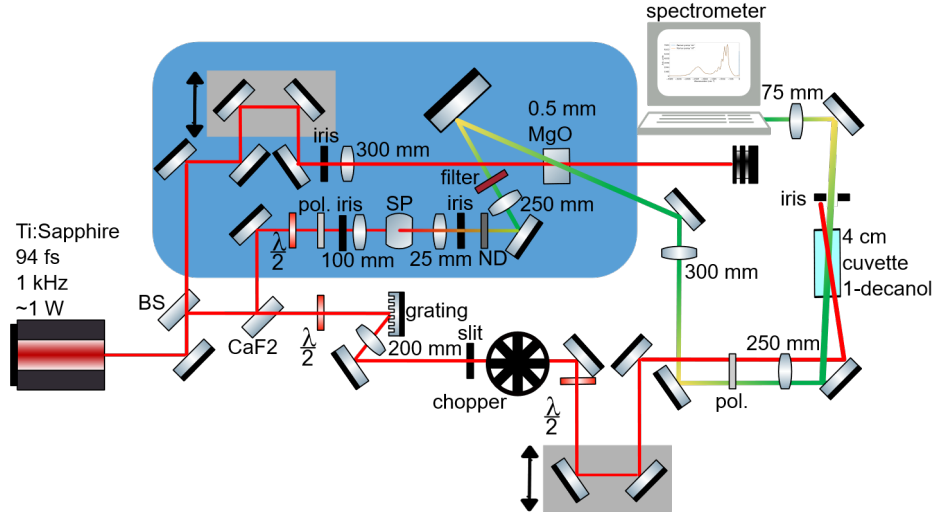


Figure 1: Experimental setup. We split the output of a Ti: Sapphire laser for Kerr-instability amplification (KIA) and generating the Raman pump pulse. The blue shaded region is the setup for KIA, where the KIA pump is overlapped temporally and spatially with the supercontinuum seed, produced by a sapphire plate (SP). We use KIA to amplify the supercontinuum and create the broadband Raman probe. The relative temporal delay of the KIA pump and seed is tuned to amplify the desired probe wavelength. The Raman pump pulse is created by limiting the laser spectrum via a grating and an adjustable slit. The Raman pump is overlapped with the Raman probe in the 4 cm sample cell which contains the 1-decanol sample. We measure the modulated Raman probe spectrum with a spectrometer.

We use a 1 W, 1 kHz Ti: Sapphire laser (Quantronix IntegraC) with a pulse duration of about 94 fs and a central wavelength of 785 nm. The beam is first split by a beamsplitter (BS). The transmitted part is the KIA pump, while the reflected part is split again by a  $\text{CaF}_2$  window. The reflection from the  $\text{CaF}_2$  is the seed for KIA, and the transmission becomes the Raman pump.

For KIA, the KIA pump passes through a delay stage to achieve temporal overlap with the seed. The power of the KIA pump is controlled to 180 mW by an iris that precedes the 300 mm lens which focuses the beam to a waist  $w_0$  of 70  $\mu\text{m}$  at the 0.5 mm MgO sample. This gives a peak intensity of approximately  $1.4 \times 10^{17} \text{ W/m}^2$  for the KIA pump, just below the damage threshold for  $\text{MgO}$ <sup>32</sup>. The MgO sample is oriented in the 110 plane to maximize amplification<sup>33</sup>.

The power of the seed for the KIA is controlled by a wire grid polarizer and halfwave plate ( $\lambda/2$ ) to be 1 mW before the 5 mm thick sapphire plate. The result is a stable supercontinuum spectrum that is amplified by the KIA pump. A 750 nm short pass filter is placed prior to the amplification process to optimize the amplification at the desired wavelength and attenuate the 785 nm part of the spectrum. The supercontinuum seed pulse is focused with a 250 mm lens and guided by a 50

mm diameter mirror to the MgO sample, where it is spatially overlapped with the KIA pump. The angle between the KIA pump and the seed is set to 5.1 degrees, enabling supercontinuum amplification<sup>28</sup>, and the temporal delay is tuned to maximize the 725 nm part of the amplified spectrum, which serves as the Raman probe.

Regarding the amount of amplification, we are operating the KIA setup below saturation, not in the maximal amplification regime. This choice is made to reduce saturation effects in the amplification process that can lead to an unstable spectral amplitude, and to limit the amount of 785 nm light that is amplified. Between the variable neutral density (ND) filter and the 750 nm cutoff filter, we achieve a stable Raman probe spectrum. While it is possible to amplify the power of the supercontinuum by 4 orders of magnitude, we choose to operate in a regime where the power of the probe at 635 nm compared to the unamplified supercontinuum is greater, on average, by two-three orders of magnitude.

The Raman pump pulse is created by spectrally limiting the original laser. The combination of a grating (Thorlabs GR25-1208, 1200 lines/mm) and an adjustable slit (Thorlabs VA100/M) significantly narrows the spectrum. The result of the reduced bandwidth is a transform-limited pulse duration increase from 94 fs to 750 fs, shown in Fig. 2 (a). The grating efficiency is optimized by rotating the polarization from vertical to horizontal polarization before the grating using a halfwave plate. With this narrower Raman pump spectrum comes improved spectral resolution, which will be discussed later alongside Fig. 5. The Raman pump is modulated by a chopper (Thorlabs MC1F10HP) that runs at 100 Hz. The laser TTL signal of 1 kHz is sent to the chopper controller (Thorlabs MC2000B), the 1 kHz signal is converted to 100 Hz to control the chopper. The 100 Hz output signal of the chopper is sent to an Arduino Uno, where it is converted to a 200 Hz signal, to account for the 500  $\mu$ s delay that it takes for the spectrometer to read an input trigger signal. The Arduino output is connected to the spectrometer (Ocean Optics Flame-S), which has an integration time of 2 ms, allowing it to measure 2 pulses “on” followed by 2 pulses “off”.

Since the Raman pump is horizontally polarized, it is rotated back to vertical polarization by a second halfwave plate. The Raman pump and probe are then passed through the same wire grid polarizer to ensure they have the same polarization. Next, the beams are passed in parallel through a 250 mm lens, focusing the pulses together within the sample. The Raman pump and Raman probe have peak intensities of approximately  $6.0 \times 10^{13} \text{ W/m}^2$  and  $4.0 \times 10^{14} \text{ W/m}^2$ , respectively, at the focus. The probe beam waist ( $w_{0\text{pr}} \approx 30 \mu\text{m}$ ) fits within the pump beam waist ( $w_{0\text{pu}} \approx 350 \mu\text{m}$ ). The peak intensity of the pump is less than that of the probe due to its longer pulse duration and larger beam size. We use a smaller probe beam waist relative to the pump to maximize our signal<sup>34</sup>. Consequently, the probe beam is more intense than the pump, but this is acceptable because the Raman loss signal scales linearly with probe intensity<sup>21,22</sup>.

If we operate our KIA setup in the maximum amplification regime of three-four orders of magnitude, we generate probe intensities on the order of  $10^{15} - 10^{16} \text{ W/m}^2$ . At such intensities we observe additional nonlinear effects in the 1-decanol sample, such as four-wave mixing and cross-phase modulation, which interfere with the Raman spectrum<sup>35</sup>. Further investigation of the optimal pump and probe intensities is needed.

The pure 1-decanol sample is contained in a glass sample cell with an optical path length of 4.0 cm, centered at the foci. The two beams diverge after the sample, and the Raman probe passes

through an iris that blocks the pump, and the probe spectrum is measured by the Ocean Optics Flame-S spectrometer. We have previously shown that the pulse duration of the probe is approximately 42 fs<sup>28</sup>, and we have previously measured 64.6 fs<sup>2</sup>/mm for the group velocity dispersion (GVD)<sup>36</sup>. Given a propagation distance of 4 cm, we calculate that the pulse duration will increase to 176 fs by propagating through the 4 cm sample cell, and this pulse duration is still much shorter than the Raman pump. We estimate the spatial overlap of the Raman pump and probe to be 2.5 cm, which is further substantiated by our observation that the signal is optimized from 2-4 cm long samples. Further discussion on pump-probe overlap will be discussed below.

Our method provides a clear advantage over spontaneous Raman spectroscopy in terms of acquisition time. We have previously measured the spontaneous Raman spectrum of 1-decanol using a 50 mW, 532 nm pen laser in a simple geometry and an integration time of 2s<sup>37</sup>. We have also measured the spontaneous Raman spectrum using a Raman microscope, where we employed a 532 nm CW laser with an average power of 32 mW, and integrated for 0.1ms 20 times. While it is possible to measure 1-decanol using either of these spontaneous Raman methods, this work is meant to illustrate a new way of creating a broadly tunable Raman probe spectrum, enabling future developments of FSRS experiments that might seek to use a variety of Raman pump wavelengths to match a resonance with a particular molecule. Although the strong  $-2900\text{ cm}^{-1}$  mode is visible within just four spectra, we acquire spectra for 2 minutes to resolve the weaker modes near  $-1500\text{ cm}^{-1}$  and improve the signal-to-noise ratio, corresponding to a total of 60,000 spectra.

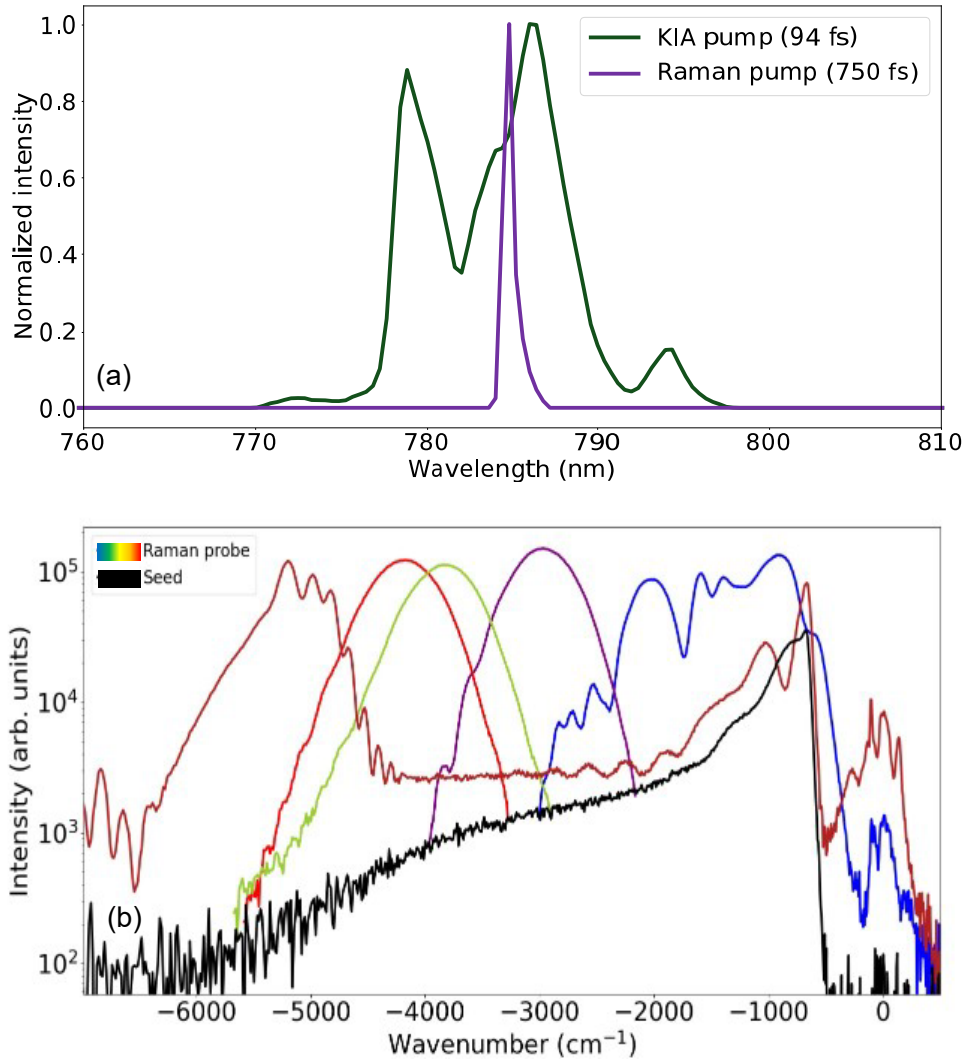


Figure 2: (a) KIA pump (green) and Raman pump (purple) spectra. The legend indicates the full width at half-maximum (FWHM) pulse durations. The laser spectral bandwidth is reduced from 9.63 nm to 1.20 nm, increasing the pulse duration from 94 fs to 750 fs. (b) adjusting the delay between the KIA pump and seed (black) allows the amplified broadband probe spectrum (various colours denoting various delays) to be tuned from  $-6000 \text{ cm}^{-1}$  to  $0 \text{ cm}^{-1}$ . The use of the 750 nm ( $-595 \text{ cm}^{-1}$ ) cutoff filter reduces amplification near  $0 \text{ cm}^{-1}$  in this case, but amplification is otherwise possible at lower energies.

The strength of our method becomes clear when considering Fig. 2 (b), which illustrates not only the level of tunability, but also the bandwidth that is available for the Raman probe. The seed spectrum is shown by the black line, and different delays are indicated by different colours in the plot. By adjusting the delay between the KIA pump and seed, the Raman probe spectrum can span from nearly  $-6000$  to  $0 \text{ cm}^{-1}$  (530 nm to 785 nm). We note once again that we are using a 750 nm ( $-595 \text{ cm}^{-1}$ ) cutoff filter prior to amplification. Without a filter, we can amplify from 533 to 1000 nm (approximately  $-6000$  to  $2800 \text{ cm}^{-1}$ ), previously not possible in a single crystal..

When a high-frequency seed pulse is amplified at around  $-5500 \text{ cm}^{-1}$  (530 nm), we observe a peak amplification of  $10^3$  and a factor of 2 amplification that extends across the supercontinuum. Given a Raman probe pulse centered at 530 nm with a bandwidth of  $2000 \text{ cm}^{-1}$ , it is possible with our method to use a short wavelength Raman pump pulse, which is advantageous for experiments that employ a pre-resonance enhancement<sup>18,38</sup>. The  $2000 \text{ cm}^{-1}$  bandwidth is not limited by the amplification process, but rather by dispersion.

## Results

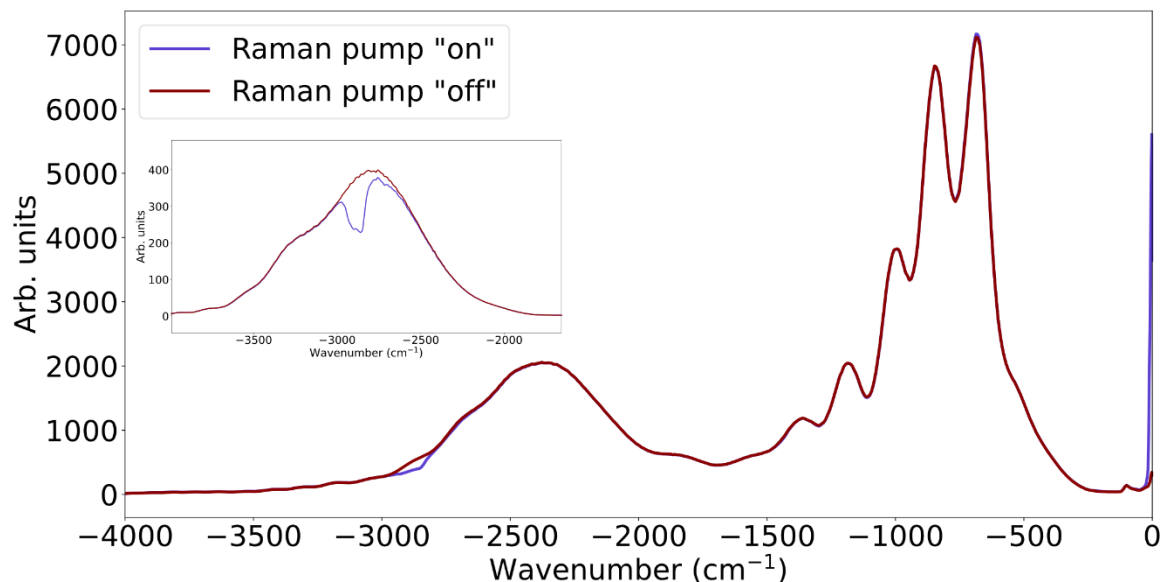


Figure 3: Raman probe spectrum with Raman pump off is shown in maroon. With the Raman pump present, stimulated inverse Raman scattering occurs, indicated by the dip in the lavender line near  $-2900 \text{ cm}^{-1}$ . The inset shows the probe spectrum when maximized at  $-2900 \text{ cm}^{-1}$ , showing the feature more clearly.

The supercontinuum (seed) spectrum resulting from sapphire is amplified by the KIA pump and optimized to span from  $-3000$  to  $-500 \text{ cm}^{-1}$  to detect the strongest Raman modes of 1-decanol. When the Raman pump pulse is included (lavender), there is a clear inverse Raman feature near  $-2900 \text{ cm}^{-1}$ .

While the stimulated Raman signal is detectable with the supercontinuum (seed) spectrum from sapphire, the chirp of the spectrum is not favourable for detecting weak signals. The chirp of a supercontinuum spectrum can be compensated for with a prism pair<sup>8</sup>. By amplifying the seed spectrum via KIA, we produce a smoother Gaussian spectrum at the desired wavenumber, improving the utility of the probe spectrum. Furthermore, amplifying the desired probe wavelength improves the signal-to-noise ratio of the measured stimulated Raman signal by a factor of 17 in our experiments.

By amplifying the supercontinuum spectrum, we generate an ultrashort, widely tunable probe spectrum. We have previously shown that it is possible to compensate for the chirp of the pulse

prior to amplification, resulting in a near transform-limited pulse<sup>31</sup>. Compressing the probe spectrum will be the subject of future work; which may facilitate detecting weaker signals.

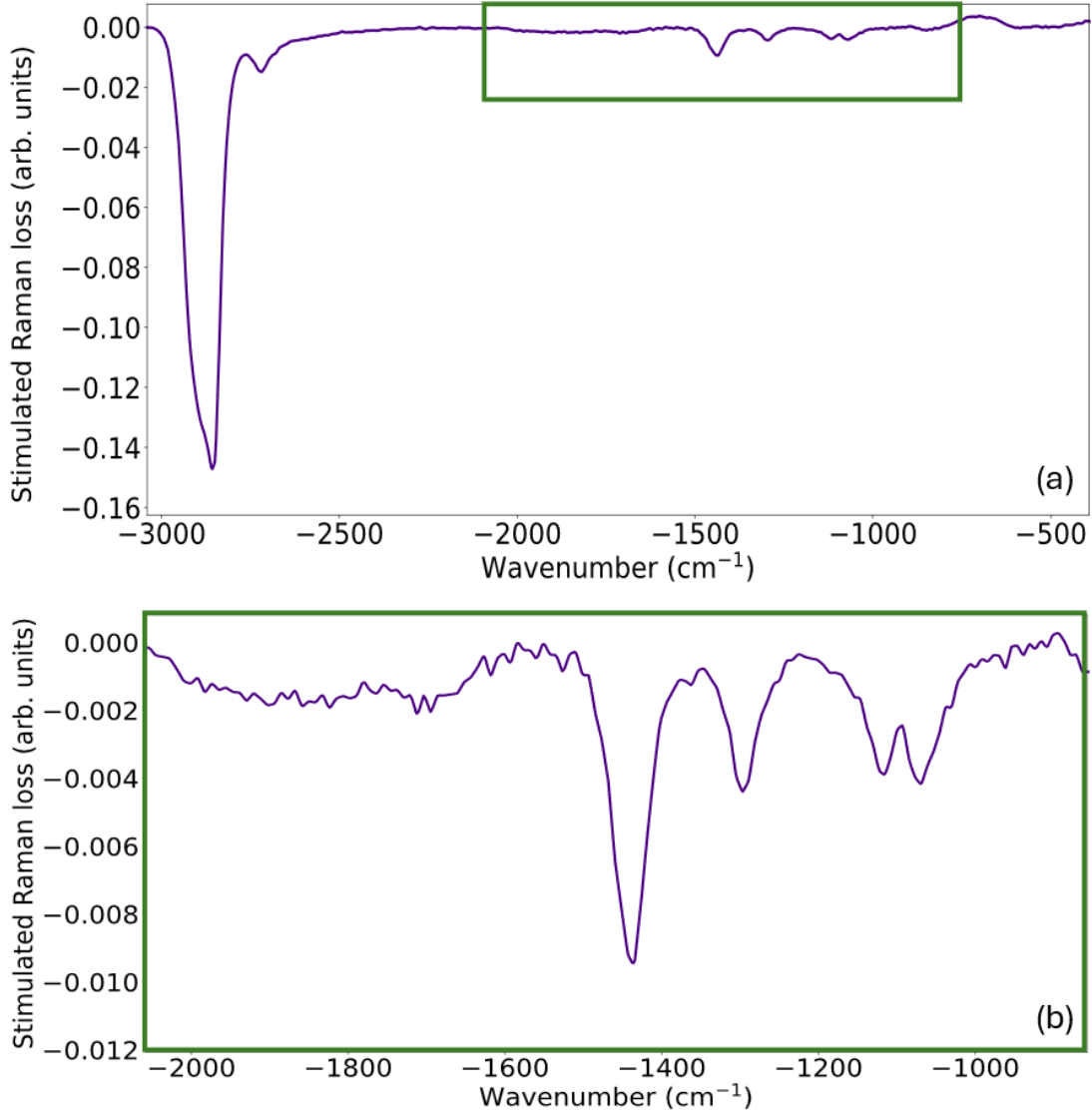


Figure 4: (a) Stimulated Raman loss spectrum calculated by dividing the blue curve by the orange curve in Fig. 3. Shown are the overtone mode at  $-2858 \text{ cm}^{-1}$ , the  $\text{CH}_2$  and  $\text{CH}_3$  stretching modes ranging from  $2880$  to  $2957 \text{ cm}^{-1}$ , and the combination mode at  $-2720 \text{ cm}^{-1}$ . (b) tuning the delay to maximize the probe spectrum at  $-1500 \text{ cm}^{-1}$ , we can resolve additional known Raman modes at lower energies.

The stimulated Raman loss spectra shown in Fig. 4 are obtained from applying the well-known equation  $\text{Raman gain (loss)} = \log_{10}\left(\frac{\text{Pump on}}{\text{Pump off}}\right)$  to the data in Fig. 3, where the measured probe



spectrum while the Raman pump is “on” is divided by the spectrum with the pump “off”<sup>2</sup>. Visually, this operation becomes obvious upon re-examining the lavender (“on”) and maroon (“off”) lines of Fig. 3. Figure 4(a) shows the full Raman loss spectrum that we measure with the probe spectrum shown previously in Fig. 3. The strongest peak near  $-2900\text{ cm}^{-1}$  is associated with the  $\text{CH}_2$  and  $\text{CH}_3$  stretching and overtone modes of 1-decanol, ranging from  $-2880$  to  $-2957\text{ cm}^{-1}$ <sup>40</sup>. The peak of the inverse Raman signal occurs at  $-2858\text{ cm}^{-1}$ , corresponding to the overtone mode<sup>40</sup>. Therefore, it is likely that we are exciting the molecule to the second excited state and driving the methyl stretching modes. In addition, the shoulder of the main peak at  $-2720\text{ cm}^{-1}$  corresponds to a known combination mode<sup>40</sup>.

Although 1-decanol has other Raman active modes at lower energies ranging from  $-1500$  to  $-1000\text{ cm}^{-1}$ , they are much weaker than those near  $-2900\text{ cm}^{-1}$ . Coupled with the chirp of our probe spectrum from  $-2000$  to  $-1000\text{ cm}^{-1}$ , characteristic of both the sapphire supercontinuum and the amount of glass used in our setup, these modes are difficult to resolve. However, by tuning the probe spectrum and maximizing the amount of light near  $-1500\text{ cm}^{-1}$ , we are in fact able to resolve these weaker modes, as shown in Fig. 4(b). More specifically, we observe the following modes:  $-1440\text{ cm}^{-1}$  ( $\text{CH}_2$  scissoring),  $-1296\text{ cm}^{-1}$  ( $\text{CH}_2$  wagging),  $-1116\text{ cm}^{-1}$  (likely  $\text{CH}_3$  rocking), and  $-1068\text{ cm}^{-1}$  which could be the CC stretching, CO stretching, or  $\text{CH}_2$  rocking mode. The OH stretching vibrations from  $-3600$  to  $-3100\text{ cm}^{-1}$  are not detected due to their reportedly low Raman cross-sections compared to the CH stretching modes<sup>40</sup>.

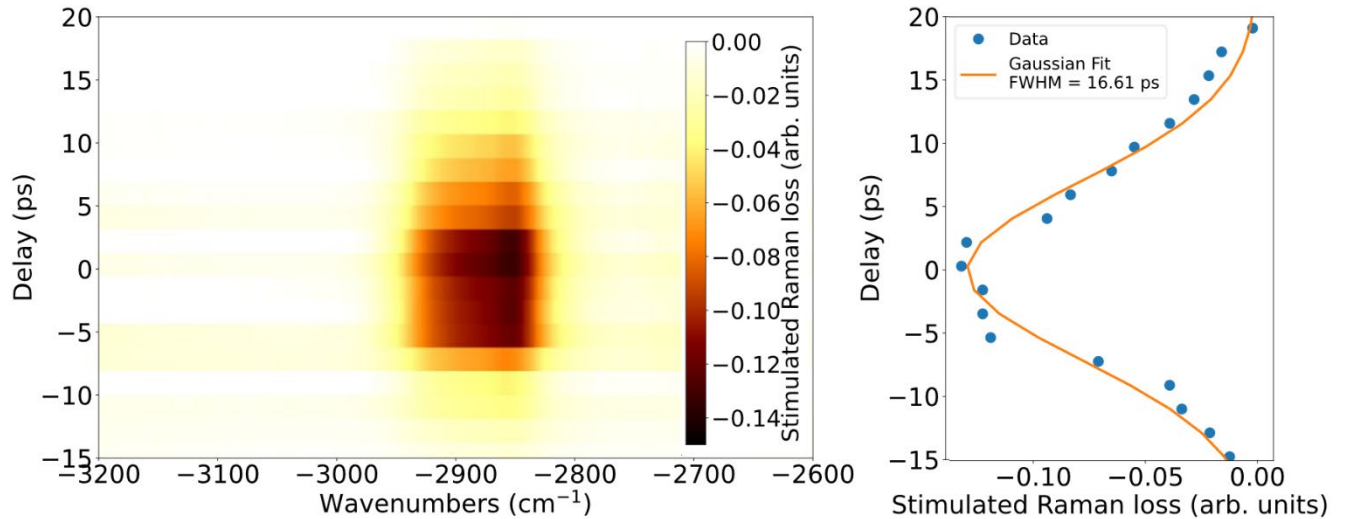


Figure 5: (a) Stimulated Raman loss as a function of pump delay (optical path difference) in picoseconds. (b) The blue points are a lineout of (a) at  $-2858\text{ cm}^{-1}$ . We fit a Gaussian function (orange line) to the data and find a full width at half-maximum (FWHM) of approximately 16.6 ps.

The spectral resolution from the FWHM pulse duration of the Raman pump is 1.2 nm, approximately equivalent to the resolution of our spectrometer. We find that the strength of the signal largely depends on the temporal overlap of the Raman pump and Raman probe pulses. Fig. 5(a) shows the Raman loss signal as a function of delay between the Raman pump and Raman probe. We convert the position of the delay stage in the Raman pump arm to optical path

difference in picoseconds. In (b), we take a vertical slice of (a) at the centre of the Raman signal, approximately  $-2858\text{ cm}^{-1}$ . We fit a Gaussian function to the data and measure a full width at half-maximum (FWHM) of approximately 16.6 ps.

## Conclusion

We have demonstrated a new method for generating a Raman probe pulse called Kerr-instability amplification for broadband stimulated Raman spectroscopy (KAB-SRS), which improves upon existing OPA systems in terms of tunability, versatility of gain media, and cost. Using KIA, we directly amplify a supercontinuum spectrum spanning from 530 to 1000 nm ( $-6000$  to  $2800\text{ cm}^{-1}$ ), enabling the use of various Raman pump wavelengths. Given our approach in having a blue-shifted probe, we measure stimulated Raman loss signals and prove our method by measuring the Raman spectrum of 1-decanol, detecting the strongest modes that occur between  $-3000$  and  $-1000\text{ cm}^{-1}$ .

With the tunability and relative affordability of our method, we expect this to be an attractive way of developing next generation FSRS setups. We anticipate our method to prove useful for detecting multiple Raman signals across a single broad spectrum.

## Acknowledgements

We thank Aaron Fisk for useful conversations. We acknowledge the technical assistance of Pratik Choudhari. We thank Steven J. Rehse for advice on the preparation of this manuscript.

## Declaration of conflicting interest

The authors declared no potential conflicts of interest with respect to the research, authorship, and/or publication of this article.

## Funding statement

N.G. Drouillard acknowledges support from the Ontario Graduate Scholarship and TJ Hammond acknowledges funding from the Natural Sciences and Engineering Research Council of Canada (RGPIN-2019-06877) and the University of Windsor Xcellerate grant (5218522). The authors also acknowledge financial support from CRC-2023-0089.

## References

1. P. Kukura, D.W. McCamant, R.A. Mathies. “Femtosecond Stimulated Raman Spectroscopy”. *Annu. Rev. Phys. Chem.* 2007. 58: 461–488.
2. D.R. Dietze, R.A. Mathies. “Femtosecond Stimulated Raman Spectroscopy”. *ChemPhysChem*. 2016. 17: 1224-1251.
3. Z. Li, M.J Deen, S. Kumar, P.R. Selvaganapathy. “Raman Spectroscopy for In-Line Water Quality Monitoring—Instrumentation and Potential”. *Sensors*. 2014. 14(9): 17275-17303.

4. M. Mulvihill, A. Tao, K. Benjauthrit, J. Arnold, P. Yang. "Surface-Enhanced Raman Spectroscopy for Trace Arsenic Detection in Contaminated Water". *Angew. Chem.* 2008. 47(34): 6456-6460.
5. D.P. Hoffman, R.A. Mathies. "Femtosecond Stimulated Raman Exposes the Role of Vibrational Coherence in Condensed-Phase Photoreactivity". *Acc. Chem. Res.* 2016. 49: 616-625.
6. C.R. Hall, J. Conyard, I.A. Heisler, G. Jones, J. Frost, W.R. Browne, B.L. Feringa, S.R. Meech. "Ultrafast Dynamics in Light-Driven Molecular Rotary Motors Probed by Femtosecond Stimulated Raman Spectroscopy". *J. Am. Chem. Soc.* 2017. 139: 7408-7414.
7. C. Fang, L. Tang, B.G. Oscar, C. Chen. "Capturing Structural Snapshots during Photochemical Reactions with Ultrafast Raman Spectroscopy: From Materials Transformation to Biosensor Responses". *J. Phys. Chem. Lett.* 2018. 9: 3253-3263.
8. D.W. McCamant, P. Kukura, R.A. Mathies. "Femtosecond Broadband Stimulated Raman: A New Approach for High-Performance Vibrational Spectroscopy". *Appl. Spectrosc.* 2003. 57(11): 1317-1323.
9. W. Peterson, J.G. de Pablo, M. Lindley, K. Hiramatsu, K. Goda. "Ultrafast impulsive Raman spectroscopy across the THz-fingerprint region". *Adv. Photonics.* 2022. 4(1): 1-12.
10. Q. Cen, Y. He, M. Xu, J. Wang, Z. Wang. "Wavelength dependent resonance Raman band intensity of broadband stimulated Raman spectroscopy of malachite green in ethanol". *J. Chem. Phys.* 2015. 142: 114201.
11. C. Fang, L. Tang. "Mapping Structural Dynamics of Proteins with Femtosecond Stimulated Raman Spectroscopy". *Annu. Rev. Phys. Chem.* 2020. 71: 239-65.
12. B.B. Oscar, C. Chen, W. Liu, L. Zhu, C. Fang. "Dynamic Raman Line Shapes on an Evolving Excited-State Landscape: Insights from Tunable Femtosecond Stimulated Raman Spectroscopy". *J. Phys. Chem. A.* 2017. 121: 5428-5441.
13. R.C. Prince, R.R. Frontiera, E.O. Potma. "Stimulated Raman Scattering: From Bulk to Nano". *Chem. Rev.* 2017. 117: 5070-5094.
14. G. Batignani, E. Pontecorvo, G. Giovannetti, C. Ferrante, G. Fumero, T. Scopigno. "Electronic resonances in broadband stimulated Raman spectroscopy". *Sci Rep.* 2016. 6: 18445.
15. P. G. Lynch, A. Das, S. Alam, C. Rich, R. R. Frontiera. "Mastering Femtosecond Stimulated Raman Spectroscopy: A Practical Guide". *ACS Phys. Chem Au.* 2024. 4 (1): 1-18.
16. C. Manzoni, G. Cerullo. "Design criteria for ultrafast optical parametric amplifiers". *J. Opt.* 2016. 18: 103501.
17. L. Zhu, W. Liu, C. Fang. "Tunable sideband laser from cascaded four-wave mixing in thin glass for ultrabroadband femtosecond stimulated Raman spectroscopy". *Appl. Phys. Lett.* 2013. 103: 061110.
18. L. Zhu, W. Liu, C. Fang. "A versatile femtosecond stimulated Raman spectroscopy setup with tunable pulses in the visible to near infrared". *Appl. Phys. Lett.* 2014. 105: 041106.
19. R.R. Frontiera, S. Shim, R.A. Mathies. "Origin of negative and dispersive features in anti-Stokes and resonance femtosecond stimulated Raman spectroscopy". *J. Chem. Phys.* 2008. 129: 064507.

20. R.E. Leighton, A.M. Alperstein, D. Punihaole, W.R. Silva, R.R. Frontiera. "Stimulated Raman versus Inverse Raman: Investigating Depletion Mechanisms for Super-Resolution Raman Microscopy". *J. Phys. Chem. B* 2023. 127: 26-36.
21. A. Aloï, R. Tommasi. "Inverse Raman Scattering in Femtosecond Broadband Transient Absorption Experiments". In: M. Khan, editor. *Raman Spectroscopy and Applications*. InTechOpen, 2017. 13. 269-290.
22. B. Mallick, A. Lakshmana, S. Umapathy. "Ultrafast Raman loss spectroscopy (URLS): instrumentation and principle". *J. Raman. Spectrosc.* 2011. 42: 1883-1890.
23. A. Barak, N. Dhiman, F. Sturm, F. Rauch, Y.A. Lakshmana, K.S. Findlay, A. Beeby, T.D. Marder, S. Umapathy. "Excited-State Intramolecular Charge-Transfer Dynamics in 4-Dimethylamino-4'-cyanodiphenylacetylene: An Ultrafast Raman Loss Spectroscopic Perspective". *ChemPhotoChem* 2022. 6: e202200146.
24. K. Roy, S. Kayal, F. Ariese, A. Beeby, S. Umapathy. "Mode specific excited state dynamics study of bis(phenylethynyl)benzene from ultrafast Raman loss spectroscopy". *J. Chem. Phys.* 2017. 146: 064303.
25. X. Qiu, X. Li, K. Niu, S. Lee. "Inverse Raman bands in ultrafast Raman loss spectroscopy". *J. Chem. Phys.* 2011. 135: 164502.
26. A. Aloï, R. Tommasi. "Inverse Raman scattering in femtosecond broadband transient absorption experiments." *Raman spectroscopy and applications*. 2017. 269-290.
27. R. Boyd. "Nonlinear Optics". Academic Press. 2020.
28. C.J. Arachchige, J.A. Stephen, TJ Hammond. "Amplification of femtosecond pulses based on  $\chi^3$  nonlinear susceptibility in MgO". *Opt. Lett.* 2021. 46 (21): 5521-5524.
29. S. Ghosh, N.G. Drouillard, TJ Hammond. "Single-stage few-cycle pulse amplification". *Phys. Rev. A*. 2024. 109: 013511.
30. S. Ghosh, N.G. Drouillard, TJ Hammond. "Supercontinuum amplification by Kerr instability". *Phys. Rev. A*. 2024. 109: 043508.
31. N.G. Drouillard, TJ Hammond. "Phase dependence of Kerr-based parametric amplification". *Phys. Rev. A*. 2024. **110**: 023517 .
32. S.Z. Xu, T.Q. Jia, H.Y. Sun, C.B. Li, X.X. Li, D.H. Feng, J.R. Qiu, Z.Z. Xu. "Mechanisms of femtosecond laser-induced breakdown and damage in MgO". *Opt. Commun.* 2006. 259: 274-280.
33. N.G. Drouillard, TJ Hammond. To be published.
34. G. Batignani, C. Ferrante, G. Fumero, M. Martinati, T. Scopigno. "Femtosecond stimulated Raman spectroscopy". *Nat. Rev. Methods Primers*. 2024. 4 (34).
35. Lim S. Lim, B. Chon, H. Rhee, M. Cho. "Spectral modulation of stimulated Raman scattering signal: Beyond weak Raman pump limit". *J. Raman Spectrosc.* 2018. 49: 607-620.
36. N. G. Drouillard, T. J. Hammond. "Measurement of dispersion and index of refraction of 1-decanol with spectrally resolved white light interferometry". *Opt. Express*. 2022. 30 (22).
37. N. G. Drouillard, "Stimulated Raman Spectroscopy with Widely Tunable Probe Pulse for Measuring Dissolved Inorganic Phosphorus in the Great Lakes," M.Sc. thesis, University of Windsor, Windsor, ON, Canada, 2024.
38. S. Laimgruber, H. Schachenmayr, B. Schmidt, W. Zinth, P. Gilch. "A femtosecond stimulated raman spectrograph for the near ultraviolet". *Appl. Phys. B*. 2006. 85: 557-564.

39. V. G. Savitski, S. Reilly, and A. J. Kemp, "Steady-State Raman Gain in Diamond as a Function of Pump Wavelength," *IEEE Journal of Quantum Electronics*, 49(2), 218–223, Feb. 2013. doi: 10.1109/JQE.2012.2237505.
40. J. Kiefer, S. Wagenfeld, D. Kerlé. "Chain length effects on the vibrational structure and molecular interactions in the liquid normal alkyl alcohols". SAA. 2018. 189: 57-65.