Compton suppressed scintillation spectrometer based on a cluster of 9×CeBr₃-NaI(Tl) phoswich detectors

M.A. Povolotskiy^{a,1}, Yu.G. Sobolev^{a,2}, S.S. Stukalov^{a,b,3}

^a Joint Institute for Nuclear Research, Dubna, Russia
 ^b Voronezh State University, Voronezh, Russia

Abstract

Measurements of the characteristics of an Compton suppressed spectrometer consisting of $9\times \text{CeBr}_3\text{-NaI}(\text{Tl})$ phoswich detectors surrounded by four CsI(Tl) scintillation detectors intended for suppressing the Compton component of the spectrum have been carried out. Measurements were performed using gamma and neutron radiation sources. Key parameters of the spectrometer have been determined: the suppression factor of the Compton part of the γ -spectrum, as well as the dependencies of the neutron detection efficiency on their energy at various detection threshold values.

PACS: 07.85.Nc, 29.40.Mc, 29.30.Hs

¹E-mail: povolotskiy@jinr.ru
 ²E-mail: sobolev@jinr.ru
 ³E-mail: stukalov@jinr.ru

Introduction

Clusters of phoswich detectors, consisting of CeBr₃-NaI(Tl) or LaBr₃-NaI(Tl) scintillators, are actively used in experiments studying the direct γ -decay processes of Giant or Pygmy Dipole Resonances (GDR, PDR) formed in reactions with charged particles [1], as well as in studies of β -decay processes of radioactive nuclei [2].

The application of phoswich detectors in modern nuclear physics experiments is due to their ability to efficiently separate high-energy γ -quanta from cascades of low-energy γ -radiation of high multiplicity.

This is especially important when analyzing events of nuclear reactions in which the total energy of the registered γ -radiation cascades is close in value to the energies of high-energy γ -quanta.

Phoswich detectors based on modern scintillators such as CeBr₃, LaBr₃, GaGG:Ce, etc. [3–5], have a high energy and time resolution. This makes them promising for applications in precision γ -spectroscopy tasks.

The high values of light output and short scintillation pulse-time of these materials, combined with the average values of the mass numbers of the nuclei composing them, also allow considering the possibility of their use as detectors for fast neutrons $E \sim 1 \div 50$ MeV.

This work is dedicated to the issues of applying a cluster consisting of $9\times \text{CeBr}_3\text{-NaI}(\text{Tl})$ phoswich detectors as an Compton suppressed γ -spectrometer, as well as a detector for fast neutrons.

The spectrometer is one of the elements of the MULTI setup [6], created at the Flerov Laboratory of Nuclear Reactions (FLNR), JINR, Dubna, and will be used both in conjunction with the 4II-spectrometer MULTI, see Fig. 1, and independently.

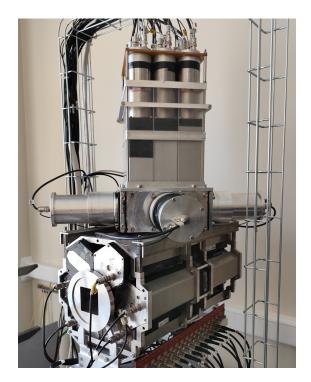


Figure 1: Photo of the MULTI experimental setup.

The spectrometer is intended for measuring total reaction cross sections with neutron-rich nuclei located far from stability line, as well as for studying their β -decay processes.

The cluster of CeBr₃-NaI(Tl) phoswich detectors as a compton suppressed spectrometer

The cluster is a close-packed assembly of nine phoswich detectors arranged in a 3×3 matrix pattern, see Fig.2b).

A phoswich detector (see Fig.2a)) consists of a pair of CeBr₃ and NaI(Tl) scintillation crystals with dimensions of $50 \times 50 \times 50 \text{ mm}^3$ and $50 \times 50 \times 150 \text{ mm}^3$, respectively, optically connected and placed in a thin-walled sealed housing with a transparent end connected to a photomultiplier.

The thickness of the CeBr₃ crystal was chosen to absorb most part of low-energy γ -quanta ($E\gamma$ less than 1 MeV) in the scintillator material, while a high-energy γ -quanta ($E\gamma$ more than 5 MeV) passed through the CeBr₃ crystal and had a first interaction with atoms of the NaI(Tl). GEANT-4 calculations showed that the probability of γ -quanta with $E\gamma = 1$ MeV passing through the CeBr₃ material into the NaI(Tl) scintillator does not exceed \sim 10%.

The CeBr₃ crystal has a single-component scintillation, characterized not only by a fast scintillation rise time (\sim 10 ns) but also a short decay time (\sim 50 ns) [7], as well as high light output, comparable in magnitude to NaI(Tl). These parameters combine excellently with the parameters of the slow NaI(Tl) scintillator within the phoswich detector to ensure effective pulse shape analysis. Fig. 3 shows a two-dimensional spectrum of the distribution of the photomultiplier anode signal charges from the fast and slow scintillation components during the registration of γ -quanta from a ⁶⁰Co source by the CeBr₃-NaI(Tl) phoswich detector.

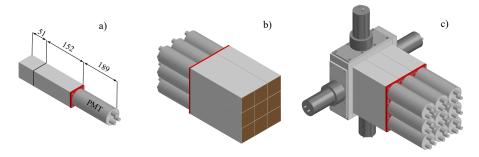


Figure 2: a) Scintillation CeBr₃-NaI(Tl) phoswich detector. b) Cluster of $9 \times$ phoswich detectors. c) Compton suppressed γ -spectrometer from a cluster of phoswich detectors and 4 suppressed Compton shield detectors.

The X-axis "FAST" represents the charge of the photomultiplier (PMT) anode pulses measured in a short time interval $\tau_{FAST} = 30$ ns. This time interval corresponds to 60% of the CeBr3 scintillation.

The Y-axis "SLOW" represents the magnitude of the pulse charge integrated in the time gate $\tau_{SLOW} = 400$ ns, which follows the "FAST"

interval. In this time gate, mainly the slow component of the NaI(Tl) scintillation was integrated.

Points populating the lower diagonal area "1" correspond to events of γ -quantum registration only in the CeBr₃ scintillator material. Points populating the area marked by contour "2" correspond to acts of γ -quantum registration only in the NaI(Tl) scintillator material. The area marked by contour "3", located between contours "1" and "2", is populated by points related to events in which the γ -quantum and its secondary interaction products, formed during the passage through the materials, were registered in both scintillators.

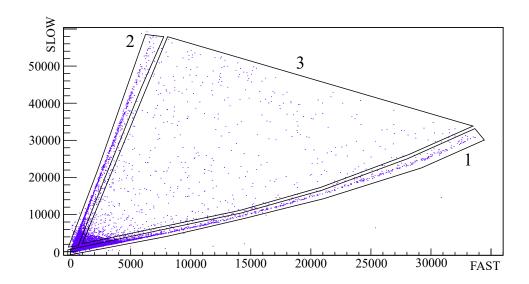


Figure 3: Two-dimensional spectrum of the distribution of PMT pulse charges, registering scintillations of CeBr₃-NaI(Tl) crystals of a phoswich detector from a ⁶⁰Co source. "1" – area of events of γ -quantum registration only in CeBr₃, "2" – area of events of γ -quantum registration only in NaI(Tl), "3" – area of events of γ -quantum registration in both CeBr₃ and NaI(Tl) scintillators

Analysis of the PMT pulse shape of the phoswich detector using the two-dimensional distribution described above allows for the separation of low-energy and high-energy γ -quanta, as well as the implementation of the Compton suppressed spectrometer mode.

We can make a conclusions about the complete absorption of the γ -quantum in the CeBr₃ crystal when two following conditions are simultaneously met, based on the analysis of γ -quantum registration events by the phoswich detector using the two-dimensional spectrum (Fig. 3) and spectra in the surrounding detectors:

-absence of signals in neighboring phoswich detectors and CsI(Tl) Compton shield detectors;

-absence of the slow scintillation component corresponding to energy deposition in the NaI(Tl) part of the phoswich detector.

Measurement of the Compton Component Suppression Factor

The study of the characteristics of the Compton suppressed spectrometer based on the cluster of phoswich detectors was carried out using γ -sources: 60 Co, 137 Cs, and 152 Eu.

Measurements of Compton suppression coefficient AC (see below) (which is a ratio of summ of the events in Compton-part of energy spectra measured in inclusive and Compton suppression modes) were performed in γ - γ coincidence using the 60 Co source. The trigger detector was a CeBr₃ scintillator. Its logical signal initiated event recording in the VME data acquisition system [8].

The Mesytec MVLC crate controller [8] received a start pulse from the trigger detector. Then the MDPP16-QDC [8] system recorded signals from all detectors of the cluster within a given time gates.

To select γ - γ coincidences, events were chosen under the conditions in which:

- the fast time pulses from the trigger and phoswich detectors are in coincidence in narrow (ΔT_{FWHM} =10ns) time window.
- the trigger CeBr₃ detector registered a total absorbed γ -quantum $E\gamma$ =1332 keV, and the cluster registered a γ -quantum with energy E_{γ} =1173 keV.

This method effectively suppresses background from the intrinsic

glow of the CeBr₃ detector.


Event selection algorithm was applied aimed at identifying the complete absorption of the γ -quantum in its CeBr₃ part of phoswich detector. The key criterion was the requirement that the γ -quantum energy be registered only in the CeBr₃ part of the phoswich detector, with a complete absence of signals in its NaI(Tl) section, as well as in neighboring detectors of the cluster and the CsI(Tl) block. This allows for the selection of events in which the Compton-scattered γ -quantum does not leave the volume of the CeBr₃ crystal.

The value of measured Peak efficiency for γ -quanta $E\gamma=1.17$ MeV (in the full absorption peak) for the CeBr₃ part of the cluster of phoswich detector, located at 10 cm from the ⁶⁰Co source, was $\varepsilon=21.8\%$. value of Peak efficiency is more than 21 times higher than one for HPGe ($\emptyset=55$ mm, H=60 mm) detectors [9] located at the same distance.

Measurements of Compton suppression coefficient with the 152 Eu γ -source were carried our in "self-trigger" mode.

¹⁵²Eu γ-source emits the cascades of γ-quanta in a wide energy range E_{γ} = 121 ÷ 1408] keV and can be used for a close approach to a real experiment on nuclear reactions, where it is impossible to unambiguously separate the Compton continuum.

Fig.4 shows examples of energy spectra of γ -sources 60 Co (a) and 152 Eu (b), measured by one of the phoswich detectors. Solid lines represent the inclusive energy spectrum. Dashed lines represent the spectrum obtained taking into account scattering into neighboring phoswich detectors. Dash-dotted lines represent the spectrum with additional consideration of scattering into the CsI(Tl) detector block. The energy threshold for Compton scattering registration was set at 60 keV.

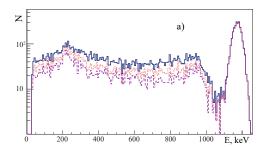

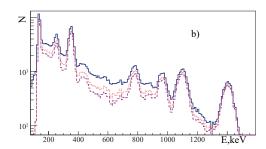


Figure 4: Energy spectra of γ -sources ^{60}Co (a) and ^{152}Eu (b) with superimposed amplitude spectra obtained under Compton suppressed shield conditions.

Data analysis with 152 Eu shows a loss of events in the full absorption peaks was observed — a known effect for Compton suppressed spectrometers. The loss of events arises due to random coincidences of two γ -quanta from one cascade, registered in neighboring detectors within the time window allotted for Compton scattering registration[10].

The reason lies in a fundamental limitation of the method: the electronic system cannot distinguish between a pair of signals from two cascade γ -quanta registered in different detectors, and a pair from one γ -quantum that underwent Compton scattering. This leads to a systematic loss of full absorption events. The severity of the effect is directly influenced by two factors: the presence of cascade transitions in the source and its intensity.

For monoenergetic sources, such as 137 Cs, or for correlated $\gamma - \gamma$ measurements with 60 Co, the shape of the spectrum allows clear separation of the energy region corresponding to Compton scattering events.

In these cases, the suppression factor AC is calculated as the ratio of summ of the events in Compton-part of energy spectra measured in inclusive and Compton suppression modes:

$$AC = (1 - \frac{N_{AC}}{N}) \cdot 100\% \tag{1}$$

where N and N_{AC} are the number of events in the Compton

region before and after applying the suppression algorithm, respectively. It is is difficult to isolate a pure Compton region of the spectrum for sources emitting cascades of γ -quanta in a wide energy range (e.g., 152 Eu) due to the overlap of numerous photo-peaks and their Compton continua. The Compton Suppression Factor (CSF) criterion "peak-to-total" [10] was applied for analyzing spectra from 152 Eu, which evaluates suppression effectiveness by the reduction of background under full absorption peaks. The calculation was performed using the following formula:

$$CSF_{total} = \left(\frac{N_{peak}}{N_{total}}\right)_{w/o.AC} / \left(\frac{N_{peak}}{N_{total}}\right)_{AC}$$
 (2)

where $\left(\frac{N_{peak}}{N_{total}}\right)_{w/o.AC}$ is the ratio of the number of events in the peak without substrate to the total number of events (total) without applying suppression, $\left(\frac{N_{peak}}{N_{total}}\right)_{AC}$ is the ratio of the number of events in the peak without substrate to the total number of events (total) with suppression applied.

Table 1 presents the AC coefficients obtained from the analysis of spectra of neighboring phoswich detectors (excluding spectra from the 4×CsI(Tl) detectors) as an example of the obtained measurement results.

Additionally, Table 1 presents the AC(CsI) suppression coefficients of the Compton part of the spectrum taking into account analysis in all detectors of the spectrometer to evaluate the effectiveness of suppressing the Compton part of the spectrum by the CsI(Tl) detectors. It is seen from Table 1 that including an additional Compton shield of $4\times \text{CsI}(\text{Tl})$ detectors increases the effectiveness of suppressing the Compton component of the spectrum by $\times 1.5$ times.

Detector	AC	AC(CsI)
	(%)	(%)
Phoswich 1	31.0 ± 1.2	47.0 ± 1.0
Phoswich 2	44.4 ± 0.9	53.3 ± 0.8
Phoswich 3	32.3 ± 1.2	45.4 ± 1.0
Phoswich 4	46.2 ± 0.9	58.0 ± 0.8
Phoswich 5	63.8 ± 0.7	65.6 ± 0.6
Phoswich 6	46.3 ± 1.0	55.1 ± 0.9
Phoswich 7	31.3 ± 1.3	50.4 ± 1.0
Phoswich 8	44.7 ± 1.0	55.8 ± 0.9
Phoswich 9	32.8 ± 1.2	49.1 ± 1.0

Table 1: Suppression factors of the Compton part of the spectrum obtained in measurements with 60 Co taking into account only the cluster of phoswich detectors (AC) and suppression by the cluster supplemented with a block of $4 \times CsI(Tl)$ detectors (AC(CsI)).

Results of measurements with subsequent processing for the γ -source 152 Eu are also presented. Due to the large volume of data because of the need to calculate suppression for the peak corresponding to each energy of γ -quanta emitted by the source, results for one detector of the cluster are presented as an example.

E_{γ}	$(N_{\rm peak}/N_{\rm tot})_{w/oAC}$	$(N_{ m peak}/N_{ m tot})_{AC}$	CSF
(keV)	$\times 10^{-3}$	$\times 10^{-3}$	
245	28.3 ± 0.3	36.3 ± 0.5	1.282 ± 0.026
344	89.9 ± 0.7	126.1 ± 1	1.402 ± 0.016
779	18.4 ± 0.3	25.6 ± 0.5	1.391 ± 0.035
964	14.7 ± 0.3	21.4 ± 0.4	1.455 ± 0.04
1112	25.8 ± 0.3	35.9 ± 0.5	1.391 ± 0.029
1408	21.1 ± 0.3	30.2 ± 0.5	1.431 ± 0.033

Table 2: Suppression coefficients of the Compton part of the spectrum obtained in measurements with the ¹⁵²Eu source with intermediate calculated peak/total values.

The comparison with modern detector systems allows highlighting the main advantages of the developed spectrometer, the basis of which is a cluster of CeBr₃-NaI(Tl) phoswich detectors with an Compton suppressed shield of four CsI(Tl) detectors.

Although the combination of an HPGe–BGO Compton suppressed clover detector remains a common solution for precision γ -spectroscopy and provides high-quality Compton background suppression, its drawbacks — low full absorption peak registration efficiency relative to CeBr₃ [9], limited HPGe detector volume, and operational complexity — make it less suitable for studying reactions with light nuclei. In such tasks, where high-energy γ -quanta prevail, the cluster Compton suppressed γ -spectrometer proves to be a more effective tool.


The measured characteristics of the phoswich detector cluster presented in this work, such as the suppression factor of the Compton component of the energy spectrum ${}_{\downarrow}AC(CsI){}_{\downarrow}\simeq 53.3\%$, energy resolution ($\Delta E(FWHM)\simeq 30~keV$), are sufficient for effectively solving a number of precision γ -spectroscopy tasks in the region of light nuclei of neutron-rich isotopes Be, B, C, N.

Neutron registration efficiency by the phoswich detector cluster

The efficiency of neutron registration by the CeBr₃ part of the phoswich detector in the energy range $E_n = 1.4\text{--}5.6$ MeV was measured using the tagged particle method. A 239 Pu/ 9 Be source was used as the neutron source, which is characterized by the fact that in 60% of neutron emission acts, the process proceeds via the 9 Be(α ,n) 12 C* channel and is accompanied by the emission of γ -quanta $E\gamma = 4.44$ MeV.

Neutrons were identified by coincidences with γ -quanta $E\gamma = 4.44$ MeV, which are emitted during the decay of the first excited state of the 12 C nucleus populated in 13 C $\rightarrow n+^{12}$ C*. Time-of-flight method(TOF) was used for identification and neutron energy determination.

Recording of an event by the VME data acquisition system was triggered by a logical signal from the trigger CeBr3 detector registering γ -quanta $E\gamma = 4.44$ MeV. Registration of a γ -quanta $E\gamma = 4.44$ MeV in the full absorption peak in the trigger detector corresponded to an act of isotropic neutron emission with a known energy distribution [11]. Events in which the trigger detector registered γ -quanta $E\gamma = 4.44$ MeV in the single and double escape peaks of annihilation quanta corresponded to the emission of both neutrons from the source and γ -quanta $E\gamma = 511$ keV from the trigger detector. Based on the known distance between the trigger detector and the phoswich detectors and the position of the γ - γ coincidence peak, the zero point of the time scale was determined, used to calculate the kinetic energy of neutrons E_n from the TOF time.

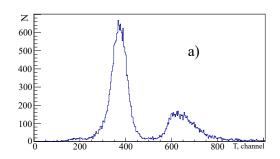

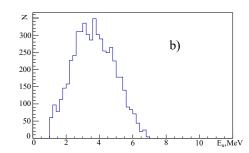


Figure 5: a) Time spectrum of γ - γ and γ – n coincidence of particles emitted from a 239 Pu/ 9 Be source; b) Energy spectrum of neutrons emitted from the 239 Pu/ 9 Be source in coincidence with γ -quanta $E\gamma = 4.44$ MeV.

The time spectrum of $\gamma - \gamma$ and $\gamma - n$ coincidences of particles emitted from the 239 Pu/ 9 Be source is shown in Fig. 5 (a). The energy spectrum of neutrons registered in coincidence with γ -quanta from the 239 Pu/ 9 Be source is shown in Fig. 5 (b). Neutron energy E_n was calculated using the non-relativistic formula $E_n = \frac{1}{2}m_nv^2$ based on time-of-flight data. The registration threshold was 60 keV on the γ -scale. The energy dependence of the neutron registration efficiency in the range $E_n = 1.4 \div 5.6$ MeV was determined by comparing the

measured neutron energy spectrum and the spectrum from the work [11], normalized to the same number of events, see Fig. 6. Fig.6 shows the energy spectra of neutrons from the 239 Pu/ 9 Be source, obtained in coincidence with γ -quanta in this work (solid histogram) and normalized to the same number of events from the work [11] (dashed histogram).

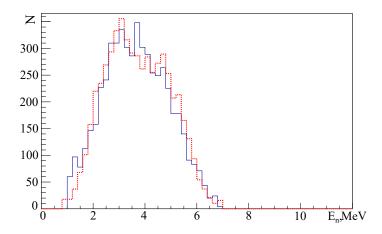
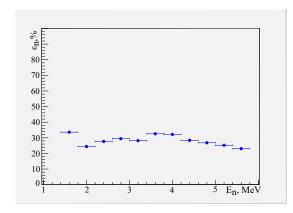


Figure 6: Energy spectrum of neutrons from this work (solid histogram). Neutron spectrum from [11] (dashed histogram).

The energy dependence of the neutron registration efficiency by the CeBr3 part of the phoswich detector cluster at different registration threshold values was obtained as a result of analyzing the neutron energy spectra. Fig.7 Left: shows the energy dependence of the neutron registration efficiency $\epsilon_n\%$ at a registration threshold Eth = 60 keV. Right: Dependence of the average registration efficiency $\bar{\epsilon}_n\%$ in the energy range $E_n = 1.4 \div 5.6$ MeV on the registration threshold value E_{th} (the abscissa axis shows the registration threshold on the energy γ -scale).



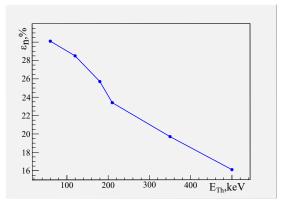


Figure 7: (Left) Energy dependence of the neutron registration efficiency $\epsilon_n(E_n)\%$ by the CeBr3 scintillator of the phoswich detector at a registration threshold Eth = 60 keV. (Right) Dependence of the average registration efficiency $j\epsilon_n\dot{\delta}\%$ on the registration threshold value E_{th} .

The highest value $\bar{\epsilon_n} \simeq 30\%$ was obtained at an energy registration threshold $E_{th} = 60$ keV. These results showed that the CeBr₃ scintillation detector allows for neutron registration with relatively high efficiency, and also confirms the main advantage of scintillators of this type — the low dependence of efficiency on the energy of the registered neutron over a wide energy range $E_n = 1.4 \div 5.6$ MeV.

Conclusion

As a result of the measurements, the main characteristics of the spectrometer based on the $9\times \text{CeBr}_3\text{-NaI}(\text{Tl})$ cluster with an Compton suppressed shield block of $4\times \text{CsI}(\text{Tl})$ were determined.

The measured characteristics of the spectrometer for γ -radiation registration (energy resolution, registration efficiency and full absorption peak efficiency, as well as suppression factors of the Compton part of the γ -quantum energy spectrum) demonstrated its high effectiveness.

Equipping the cluster with four CsI(Tl) detectors allowed increasing the suppression factor of the Compton part of the spectrum by

a maximum of ~ 1.5 times.

The maximum suppression factor value was $\sim 65\%$ for gamma radiation from the $^{60}\mathrm{Co}$ source.

The dependence of registration efficiency on energy $\varepsilon(En)$ in the range $En=1.4\div5.6$ MeV when registering neutrons from the ²³⁹Pu/⁹Be source was obtained. The efficiency demonstrates a flat dependence on energy with an average value $\varepsilon(E_n)=\simeq30\%$ for the entire range, which agrees with data previously obtained for inorganic scintillators [12].

Additionally, the dependence of efficiency on the energy registration threshold was investigated, demonstrating the necessity of operating at minimal thresholds. The authors thank the Russian Science Foundation for financial support (Project N 24-22-00117).

References

- 1. γ decay to the ground state from the excitations above the neutron threshold in the Pb 208 (p, p' γ) reaction at 85 MeV / B. Wasilewska [et al.] // Physical Review C. 2022. Vol. 105, no. 1. P. 014310.
- 2. β -delayed spectroscopy of Ge 48 80 and competition between Gamow-Teller and first-forbidden transitions in Ga 49 80 g+ m β decay / R. Li [et al.] // Physical Review C. 2025. Vol. 112, no. 4. P. 044306.
- 3. Moses W. W., Shah K. S. Potential for RbGd2Br7: Ce, LaBr3: Ce, LaBr3: Ce, LaBr3: Ce in nuclear medical imaging // Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment. 2005. Vol. 537, no. 1/2. P. 317–320.
- 4. First performance results of Ce: GAGG scintillation crystals with silicon photomultipliers / J. Y. Yeom [et al.] // IEEE transactions on nuclear science. 2013. Vol. 60, no. 2. P. 988–992.

- 5. Scintillation and detection characteristics of high-sensitivity CeBr₃ γ -ray spectrometers / F. Quarati [et al.] // Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment. 2013. Vol. 729. P. 596–604.
- 6. MULTI-2, a 4π spectrometer for total reaction cross section measurements / I. Siváček [et al.] // Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment. — 2020. — Vol. 976. — P. 164255.
- 7. Characteristics and time resolutions of two CeBr3 gamma-ray spectrometers / J. Qin [et al.] // Radiation Detection Technology and Methods. 2020. Vol. 4, no. 3. P. 327–336.
- 8. https://www.mesytec.com/.
- 9. Challan M. B. Gamma-ray efficiency of a HPGe detector as a function of energy and geometry // Applied Radiation and Isotopes. 2013. Vol. 82. P. 166–169.
- 10. Monte-Carlo optimisation of a Compton suppression system for use with a broad-energy HPGe detector / R. Britton [et al.] // Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment. 2014. Vol. 762. P. 42–53.
- 11. A comparison of untagged gamma-ray and tagged-neutron yields from 241 AmBe and 238 PuBe sources / J. Scherzinger [et al.] // Applied Radiation and Isotopes. 2017. Vol. 127. P. 98–102.
- 12. Response of BaF₂, CsI (Tl) and Pb-glass detectors to neutrons below 22 MeV / T. Matulewicz [et al.] // Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment. 1989. Vol. 274, no. 3. P. 501–506.